


Inspired by temperature

Betriebsanleitung · Operation manual · Manual de instrucciones · Manuel d'utilisation · Manuale de d'uso · **사용 설명서** · Manual de instruções · Инструкция по эксплуатации · Kullanım talimatı · 操作说明书 ·

Betriebsan instruccion · **사용 설명** по эксплуа Betriebsan instruccion

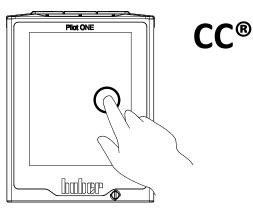

Manual de ale de d'uso нструкция **: 说明书** Manual de ale de d'uso


· 사용 설명서 · Manual de instruções · Инструкция по эксплуатации · Kullanım talimatı · 操作说明书

This documentation does not contain a device-specific technical appendix.

You can request the full operating instructions from info@huber-online.com. Please give the model designation and serial number of your temperature control unit in your e-mail.





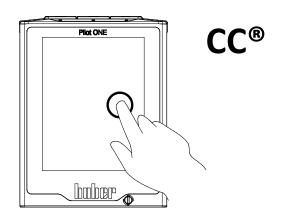



CC®

Pilot ONE®

This operation manual is a translation of the original operation manual.


VALID FOR:

K6 K1x K2x CC®-K6 CC®-Kxx Variostat® Ministat® 125 Ministat® 2xx CC®-4xx CC®-5xx CC®-5xx CC®-8xx CC®-8xx

Abbreviations used in model names: without = air cooled, S = higher cooling capacity, w = water cooled, wl = water/air cooled

5



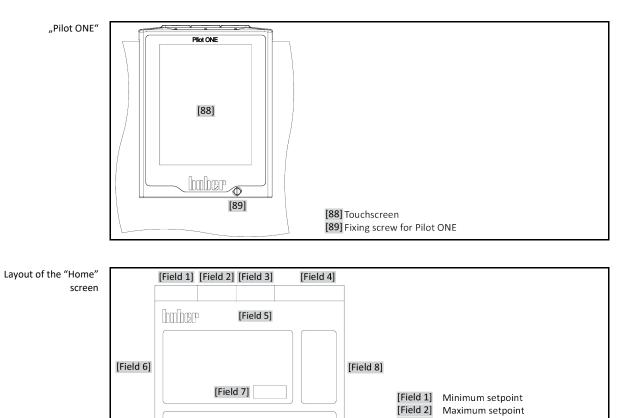

[Field 3] Set over-temperature protection

[Field 7]Button for virtual keyboard[Field 8]Level and pump status[Field 9]Process graphic

[Field 5] Date and time [Field 6] Process, internal and setpoint value

[Field 4] Status pictogram

[Field 10] [Field 10] Status cell [Field 11] [Field 11] Touchbuttons

[Field 9]

CC®

6

Table of contents

V2.5.0en/27.07.21//17.12

1	Introduction 14
1.1	Identification / symbols in the operation manual14
1.2	Information on the EU Declaration of Conformity14
1.3	Safety
1.3.1	Symbols used for Safety Instructions14
1.3.2	Representation of safety identifiers on the temperature control unit
1.3.3	Proper operation
1.3.4	Reasonably foreseeable misuse16
1.4	Responsible bodies and operators – Obligations and requirements
1.4.1	Obligations of the responsible body
1.4.1.1	Proper disposal of resources and consumables
1.4.1.2	Temperature control unit with natural refrigerants (NR)
1.4.1.3	Temperature control units with fluorinated greenhouse gases/refrigerants20
1.4.2	Requirements for operators
1.4.3	Obligations of the operators
1.5	General information
1.5.1	Description of workstation20
1.5.2	Safety devices to DIN 1287621
1.5.2.1	Temperature control unit (without heater) with classic float switch21
1.5.2.2	Temperature control units (with heater) with electronic low level sensor21
1.5.3	Further protective devices
1.5.3.1	Power interruption
1.5.3.2	Alarm functions
1.5.3.3	Warning messages22
1.6	Exemplary illustrations of the cooling variants
1.6.1	Consequence of inadequate energy dissipation
2	Commissioning
2	Commissioning 25
2.1	In-plant transport
2.1.1	Lifting and transporting the temperature control unit
2.1.1.1	Temperature control unit with lifting eyes
2.1.1.2	Temperature control unit without lifting eyes
2.1.2	Mounting/removing leveling feet
2.1.3	Positioning the temperature control unit27
2.1.3.1	Temperature control unit with casters27
2.1.3.2	Temperature control unit without casters
2.2	Unpacking
2.3	Ambient conditions
2.3.1	EMC-specific notes
2.4	Installation conditions
2.5	Recommended temperature control and cooling water hoses
2.6	Wrench sizes and torques
2.7	Temperature control units with water cooling
2.8	Temperature control unit with batch
2.8.1	Operation as bath thermostat
2.9	Preparations for operation
2.9.1	Unscrewing/activating the leveling feet (if any)



2.9.2	Opening/closing valves	34
2.9.3	Installing collecting container	34
2.9.4	Cooling/heating baths: Connect the energizing lead	34
2.9.5	Connecting the functional earth	35
2.10	Connecting externally closed application	35
2.10.1	Connecting an externally closed application	35
2.11	Connecting to the power supply	35
2.11.1	Connection using socket with protective earth (PE)	35
2.11.2	Connection via hard wiring	36
2.11.3	Converting the power supply connection	36
2.11.4	Cooling/heating baths: Connecting to the power supply	37
2.11.4.1	Utilization as a pure cooling bath (without hanger thermostat)	37
2.11.4.2	Utilization as a cooling/heating baths (with hanger thermostat)	37
3	Function description	39
3.1	Function description of the temperature control unit	39
3.1.1	General functions	39
3.1.2	Other functions	39
3.2	Information on the thermal fluids	40
3.2.1	Only valid for CC [®] -520w, CC [®] -525w, CC [®] -906w and cooling baths used i	
	continuous operation	
3.3	To be noted when planning the test	
3.4	"Pilot ONE®" controller	
3.4.1	Functional overview of "Pilot ONE®"	
3.5	Clock/event function	
3.5.1	Rechargeable accumulator	
3.5.2	Programmable event function	
3.5.2.1	Event function "Alarm clock event"	
3.5.2.2	Event function "Program event"	
3.6	Operation via the touch screen	
3.7	Display instruments	
3.7.1	The touchscreen [88]	45
3.8	Control instruments	46
3.8.1	The touchbuttons	46
3.8.2	The categories	46
3.8.3	The sub-categories	46
3.8.4	The dialogs	46
3.9	Function examples	
3.9.1	Display of software version	
3.9.2	Start & Stop	
3.9.3	Copying the settings to a data carrier	
3.9.3.1	Saving to a USB flash drive	48
3.9.3.2	Loading from a USB flash drive	48
3.9.4	Restore factory settings	
3.9.4.1	Restore to factory settings without overtemperature protection	50
3.9.4.2	Restore to factory settings including overtemperature protection	50
4	Setup mode	51
4.1	Setup mode	51
4.1.1	Turning on the temperature control unit	51
4.1.2	Turning off the temperature control unit	51
4.1.3	Cooling/heating baths: Power On/Off	52



4.1.3.1	Cooling/heating baths: Switching on (without an immersion circulator thermostat)	52
4.1.3.2	Cooling/heating baths: Switching off (without an immersion circulator thermostat)	
4.1.3.3	Cooling/heating baths: Switching on (with an immersion circulator thermostat)	52
4.1.3.4	Cooling/heating baths: Switching off (with an immersion circulator thermostat)	53
4.1.4	Setting the overtemperature protection	53
4.1.4.1	General information on the overtemperature protection	53
4.1.4.2	Setting "OT limit: heating"	54
4.1.4.3	Setting "OT expansion vessel"	
4.1.4.4	Setting "Process Safety"	
4.1.4.5	Monitoring via "Display OT values"	55
4.1.5	Testing overtemperature protection for functionality	55
4.1.6	Adjusting the Delta T limiter	55
4.1.6.1	Changing the Delta T limiter	56
4.2	The temperature control circuit	56
4.2.1	Select temperature control: Internal or process	56
4.2.2	Temperature control to internal temperature	57
4.2.3	Temperature control to process temperature	57
4.2.4	Delta T limiter	58
4.2.5	Monitoring the Pt100 temperature sensors	58
4.2.6	Optimum control parameters for optimum temperature control	58
4.2.7	Sub-category: "Select auto/expert mode"	
4.2.8	Sub-category: "Configuration auto"	59
4.2.8.1	Sub-category: "Find parameters"	59
4.2.8.2	Sub-category: "Control Dynamics"	60
4.2.8.3	Sub-category: "Fluid Properties"	62
4.2.8.4	Sub-category: "Display parameters"	63
4.2.9	Sub-category: "Configuration manual"	
4.2.9.1	Sub-category: "Change parameters"	63
4.2.9.2	Sub-category: "Display parameters"	64
4.2.9.3	Sub-category: "Control structure"	
4.2.10	Sub-category: "Reset parameters"	
4.2.11	Sub-category: "Display parameters"	65
4.2.12	Setting the setpoint thresholds	65
4.2.13	Setting the setpoint	
4.3	Filling, venting, degassing and draining	
4.3.1	Filling, venting, degassing and draining the bath thermostat	67
4.3.1.1	Filling and venting the bath thermostat and the externally closed application	67
4.3.1.2	Degassing of bath thermostat	
4.3.1.3	Draining the bath thermostat	69
5	Normal operation	70
5.1	Automatic operation	70
5.1.1	Temperature control	70
5.1.1.1	Starting the temperature control process	70
5.1.1.2	Ending the temperature control process	70
5.1.2	Temperature control via a created temperature control program	71
5.1.2.1	Starting the temperature control program	
5.1.2.2	Ending/cancelling the temperature control program	71

cc® 9

huber

OPERATION MANUAL

5.2	Cooling/heating baths: Automatic operation (without an immersion circ	ulator
F 2 4	thermostat)	
5.2.1 5.2.1.1	Cooling/heating baths: Temperature control	
-	Cooling/heating baths: Starting the temperature control process	
5.2.1.2	Cooling/heating baths: Ending the temperature control process	72
6	Interfaces and software update	73
6.1	Interfaces at the "Pilot ONE®" controller	73
6.1.1	10/100 Mbps Ethernet for RJ45 network sockets	73
6.1.2	USB-2.0 interface	73
6.1.2.1	USB-2.0 interface, host	74
6.1.2.2	USB-2.0 interface, device	
6.2	Cooling/heating baths: Interfaces at the back	74
6.2.1	Activation connector	
6.3	Cooling/heating baths: Interfaces on the hanger thermostat	
6.3.1	Service interface	
6.3.2	RS232 (option) / activation of cooling bath	
6.3.3	Connection jack for Pt100 process control sensor	
6.4	Interfaces at the Com.G@te [®] (optional)	
6.4.1	Jack LEVEL (Com.G@te [®] external only)	
6.4.2	Connector POKO (floating contact) alarm	
6.4.3	Jack AIF Reg-E-Prog	
6.4.4	Jack ECS (External Control Signal) standby	
6.4.5	Jack RS232/RS485 serial	79
6.5	Firmware update	79
7	Service/maintenance	80
7.1	Messages from the temperature control unit	80
7.1 7.2	Messages from the temperature control unit Electrical fuse	
	Electrical fuse	80
7.2		80 80
7.2 7.3	Electrical fuse Replacing the "Pilot ONE®"	80 80 81
7.2 7.3 7.4	Electrical fuse	80 80 81 81
7.2 7.3 7.4 7.4.1	Electrical fuse Replacing the "Pilot ONE®" Maintenance Function check and visual inspection	80 80 81 81 82
7.2 7.3 7.4 7.4.1 7.4.2	Electrical fuse Replacing the "Pilot ONE®" Maintenance Function check and visual inspection Replacing temperature control or coolant hoses	80 80 81 81 81 82 82 82
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1	Electrical fuse	80 80 81 81 81 82 82 83
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2	Electrical fuse	80 80 81 81 82 82 82 83 83
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3	Electrical fuse	80 80 81 81 82 82 83 83 83 84
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4	Electrical fuse	80 80 81 81 82 82 83 83 nit).84 85
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5	Electrical fuse	80 81 81 82 82 83 83 nit).84 85
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1	Electrical fuse	80 81 81 82 82 83 83 83 83 85 85
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2	Electrical fuse	80 80 81 81 82 82 83 83 83 85 85 87
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6	Electrical fuse	80 80 81 81 82 82 83 83 83 85 85 85 87 87
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6 7.7	Electrical fuse	80 81 81 82 82 83 83 83 85 85 85 87 87 88
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.1 7.5.2 7.6 7.7 7.8	Electrical fuse	80 81 81 82 82 83 83 83 85 85 85 87 87 88
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6 7.7 7.8 7.9	Electrical fuse	80 80 81 81 82 82 83 83 83 85 85 85 85 87 87 88 88 88 88 88
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6 7.7 7.8 7.9 8	Electrical fuse	80 81 81 82 82 83 83 83 85 85 85 87 87 87 88 88 88 88 89 89
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6 7.7 7.8 7.9 8 8.1	Electrical fuse Replacing the "Pilot ONE®" Maintenance Function check and visual inspection Function check and visual inspection Replacing temperature control or coolant hoses Replacing temperature control hoses Replacing coolant hoses Clean liquefier fins (air-cooled temperature control unit) Clean hat-type strainer (dirt trap) (water-cooled temperature control unit) Clean hat-type strainer (dirt trap) (water-cooled temperature control unit) Thermal fluid inspection, replacement and circuit cleaning Thermal fluid inspection Rinsing the thermal fluid circuit Cleaning the surfaces Inspect the mechanical seal Plug contacts Decontamination/repairs Shutting down Safety instructions and basic principles	80 81 81 82 82 83 83 85 85 85 87 87 87 88 88 88 88 89 89 89
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6 7.7 7.8 7.9 8 8.1 8.2	Electrical fuse	80 80 81 81 82 82 83 83 83 85 85 85 87 87 87 88 88 88 89 89 90
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2 7.4.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6 7.7 7.8 7.9 8 8.1 8.2 8.3	Electrical fuse Replacing the "Pilot ONE®" Maintenance Function check and visual inspection Replacing temperature control or coolant hoses Replacing temperature control hoses Replacing coolant hoses Clean liquefier fins (air-cooled temperature control unit) Clean hat-type strainer (dirt trap) (water-cooled temperature control unit) Clean hat-type strainer (dirt trap) (water-cooled temperature control unit) Thermal fluid inspection, replacement and circuit cleaning Thermal fluid inspection Rinsing the thermal fluid circuit Cleaning the surfaces Inspect the mechanical seal Plug contacts Decontamination/repairs Shutting down Safety instructions and basic principles Switch-off Draining the temperature control unit Draining the cooling water	80 80 81 81 82 82 83 83 83 85 85 85 85 87 87 88 88 88 89 89 90 90
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2 7.4.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6 7.7 7.8 7.9 8 8.1 8.2 8.3 8.4	Electrical fuse Replacing the "Pilot ONE®" Maintenance Function check and visual inspection Replacing temperature control or coolant hoses Replacing temperature control hoses Replacing coolant hoses Clean liquefier fins (air-cooled temperature control unit) Clean hat-type strainer (dirt trap) (water-cooled temperature control u Thermal fluid inspection, replacement and circuit cleaning Thermal fluid inspection Rinsing the thermal fluid circuit Cleaning the surfaces Inspect the mechanical seal Plug contacts Decontamination/repairs Shutting down Safety instructions and basic principles Switch-off Draining the temperature control unit	80 81 81 82 82 83 83 83 85 85 85 85 87 87 87 87 88 88 88 89 90 90
7.2 7.3 7.4 7.4.1 7.4.2 7.4.2.1 7.4.2.2 7.4.3 7.4.4 7.5 7.5.1 7.5.2 7.6 7.7 7.8 7.9 8 8.1 8.2 8.3 8.4 8.4.1	Electrical fuse	80 81 81 82 82 83 83 85 85 85 85 87 87 87 88 88 88 89 90 90 90



8.7	Packing	
8.8	Shipping	
8.9	Disposal	
8.10	Contact data	
8.10.1	Telephone number: Customer Support	92
8.10.2	Telephone number: Sales	92
8.10.3	Email address: Customer Support	92
8.11	Certificate of Compliance	93
9	Annex	94





Foreword

Dear Customer,

Thank you for choosing a temperature control unit from Peter Huber Kältemaschinenbau SE. You have made a good choice. Thank you for your trust.

Please read the operation manual carefully before putting the unit into operation. Strictly follow all notes and safety instructions.

Follow the operation manual with regard to transport, start-up, operation, maintenance, repair, storage and disposal of the temperature control unit.

We fully warrant the temperature control unit for the specified intended operation.

The models listed on page 5 are referred to in this operation manual as temperature control units and Peter Huber Kältemaschinenbau SE as Huber company or Huber.

Liability for errors and misprints excluded.

The following trademarks and the Huber logo are registered trademarks of Peter Huber Kältemaschinenbau SE in Germany and/or other countries worldwide: BFT®, CC®, Chili®, Com.G@te®, Compatible Control®, CoolNet®, DC®, E-grade®, Grande Fleur®, Huber Piccolo®, KISS®, Minichiller®, Ministat®, MP®, MPC®, Peter Huber Minichiller®, Petite Fleur®, Pilot ONE®, RotaCool®, Rotostat®, SpyControl®, SpyLight®, Tango®, TC®, UC®, Unical®, Unichiller*, Unimotive*, Unipump*, Unistat*, Unistat Tango*, Variostat*. The following trademarks are registered in Germany to DWS Synthesetechnik: DW-Therm®, DW-Therm HT®. The following trademark is a registered trademark of BASF SE: Glysantin[®].

Introduction 1

1.1 Identification / symbols in the operation manual

The following identifications and symbols are used in the texts and illustrations.

Overview	Identification / symbol Description			
	\rightarrow	Reference to information / procedure.		
	»TEXT«	Reference to a chapter in the operation manual. In the digital version, the text is clickable.		
	>TEXT< [NUMBER]	Reference to the wiring diagram in the annex. The designation and the search digit are specified.		
	>TEXT< [LETTER]	Reference to a drawing in the same paragraph. The designation and the search digit are specified.		
	•	List, first level		
	-	List, second level		

Information on the EU Declaration of Conformity 1.2

The equipment complies with the basic health and safety requirements of the European Directives CE listed below:


- Machinery Directive
- Low Voltage Directive
- EMC Directive

1.3 Safety

1.3.1 Symbols used for Safety Instructions

Safety instructions are marked by the below combinations of pictograms and signal words. The signal word describes the classification of the residual risk when disregarding the operation manual.

DANGER	Denotes an immediate hazardous situation that will result in death or serious injuries.
WARNING	Denotes a general hazardous situation that may result in death or serious injuries.
	Denotes a hazardous situation that can result in injury.
ΝΟΤΕ	Denotes a situation that can result in property material damage.
INFORMATION	Denotes important notes and usable hints.

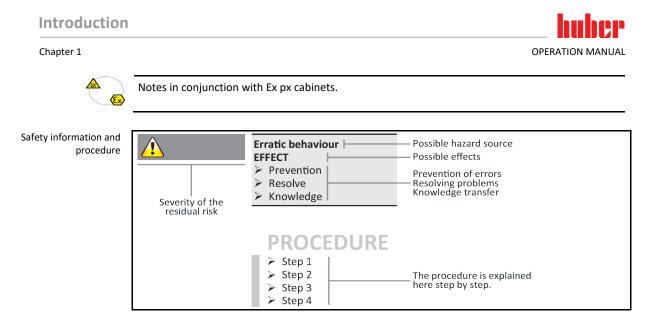

The safety information in this operation manual is designed to protect the operating company, the operator and the equipment from damage. First inform yourself about any residual risks due to misuse before you start an operation.

1.3.2 Representation of safety identifiers on the temperature control unit

The following pictograms are used as safety identifiers. The table gives an overview of the safety identifiers used here.

Overview	Identifier	Description			
	Mandatory sign				
		- Observe the instructions			
	Warning sign				
		 General warning sign Observe the instructions 			
		- Warning of electrical voltage			
		- Warning of hot surface			
		- Warning of flammable substances			

1.3.3 Proper operation

Operating the temperature control unit in a potentially explosive area DEATH THROUGH EXPLOSION

 > Do NOT install or start up the temperature control unit within an ATEX zone.

15



Chapter 1

	Improper use
	SERIOUS INJURY AND PROPERTY DAMAGE
	 Store the operation manual where it is easy to access in close proximity to the temperature control unit. Only adequately qualified operators may work with the temperature control unit. Operators must be trained before handling the temperature control unit. Check that the operators have read and understood the operation manual. Define precise responsibilities of the operators. Personal protective equipment must be provided to the operators.
	Be sure to follow the responsible body's safety rules to protect life and limb and to limit dam ages!
ΟΤΕ	
ΟΤΕ	ages!

The temperature control unit is made for industrial use. The temperature control unit is used to maintain the temperature of applications, such as glass or metal reactors or other expedient items in laboratories and industry. Flow-through coolers and calibration baths must be used only in combination with Huber temperature control units. Only use thermal fluids suitable for the overall system. The cooling or heating capacity is provided at the pump connections or - where present - in the tempering bath. For the technical specification, refer to the datasheet. \rightarrow From page 94, section **»Annex«**. Install, set up and operate the temperature control unit according to the instructions in this operation manual. Any failure to comply with the operation manual is considered as improper operation. The temperature control unit was manufactured according to the state of the art and the recognized safety rules and regulations. Safety devices are installed in your temperature control unit.

1.3.4 Re

Reasonably foreseeable misuse

Without an Ex px cabinet, the temperature control unit / accessory is **NOT** protected against explosion and must **NOT** be installed or put into operation within an ATEX Zone. When operating the temperature control unit /accessory in conjunction with an Ex px cabinet, the information in the annex (Section ATEX operation) must be observed and followed. This annex is only provided for temperature control units /accessories delivered with an Ex px cabinet. If this annex is missing, please immediately contact the Customer Support. \rightarrow Page 92, section **»Contact data«**.

Use with medical devices (e.g. in Vitro diagnostic procedure) or for direct foodstuff temperature control is **NOT** permissible.

The temperature control unit must **NOT** be used for any purposes other than temperature control in accordance with the operation manual.

The manufacturer accepts **NO** liability for damage caused by **technical modifications** to the temperature control unit, **improper handling** or use of the temperature control unit if the operation manual is **not observed**.



1.4 Responsible bodies and operators – Obligations and requirements

1.4.1 Obligations of the responsible body

The operation manual is to be stored where it is easy to access in close proximity to the temperature control unit. Only adequately qualified operators (e.g. chemists, CTA, physicists etc.) are permitted to work with the temperature control unit. Operators must be trained before handling the temperature control unit. Check that the operators have read and understood the operation manual. Define precise responsibilities of the operators. Personal protective equipment must be provided to the operators.

- The responsible body must install a condensation water / thermal fluid drip tray below the temperature control unit.
- The use of a drip tray may be prescribed by national legislation for the installation area of the temperature control unit (incl. accessory). The responsible body must check and apply the national regulations applicable for it accordingly.
- The temperature control unit complies with all applicable safety standards.
- Your system, which uses our temperature control unit, must be equally safe.
- The responsible body must design the system to ensure it is safe.
- Huber is not responsible for the safety of your system. The responsible body is responsible for the safety of the system.
- Although the temperature control unit provided by Huber meets all the applicable safety standards, integration into a system may give rise to hazards that are characteristic of the other system's design and beyond the control of Huber.
- It is the responsibility of the system integrator to ensure that the overall system, into which this temperature control unit is integrated, is safe.
- The >Mains isolator< [36] (if present) can be locked in the off position to facilitate safe system
 installation and maintenance of the temperature control unit. It is the responsibility of the responsible body to develop any lock-out/tag-out procedure for the energy source in accordance
 with local regulations (e.g. CFR 1910.147 for the US).

1.4.1.1 Proper disposal of resources and consumables

Do comply with all national disposal regulations applicable for you. Contact your local waste management company for any questions concerning disposal.

Overview	Material / Aids	Disposal / Cleaning
	Packaging material	Keep the packaging material for future use (e.g. transport).
	Thermal fluid	Please refer to the safety data sheet of the thermal fluid used for information on its proper disposal. Use the original thermal fluid container when disposing it.
	Filling accessories, e.g. beaker	Clean the filling accessories for reuse. Make sure that the materials and cleaning agents used are properly disposed of.
	Aids such as towels, cleaning cloths	Tools used to take up spilled thermal fluid must be disposed of in the same fashion as the thermal fluid itself. Tools used for cleaning must be disposed of depending on the cleaning agent used.
	Cleaning agents such as stainless steel cleaning agents, sensitive-fabrics detergents	Please refer to the safety data sheet of the cleaning agent used for information on its proper disposal. Use the original containers when disposing of large quantities of cleaning agents.
	Consumables such as air filter mats, temperature control hoses	Please refer to the safety data sheet of the consumables used for information on their proper disposal.



Chapter 1

1.4.1.2 Temperature control unit with natural refrigerants (NR)

Over 8 g refrigerant per m³ room air DEATH OR SERIOUS INJURY DUE TO EXPLOSION

- Observe the rating plate (amount of natural refrigerant contained) and the room size (maximum room concentration of natural refrigerant in case of leakage) when installing the temperature control unit.
- Over 8 g refrigerant per m³ room air: A gas warning sensor must be fitted and functioning.
- The gas warning sensor must be calibrated and maintained at regular intervals (between 6 and 12 months)
- The temperature control unit is not approved for operation in an ATEX zone.

Huber products with natural refrigerants work with numerous proven, safe and highly-sustainable technologies. The relevant standards and regulations for temperature control units with natural refrigerants contain a number of stipulations, the importance of complying with which is set out below. Please additionally: \rightarrow Page 15, section **»Proper operation**«.

Huber temperature control units are constructed to be permanently sealed and are carefully checked for leak tightness. Temperature control units with more than 150 g natural refrigerant are equipped with an additional gas warning sensor. To find out whether your temperature control unit is equipped with a gas warning sensor, refer to the data sheet. \rightarrow From page 94, section **»Annex**«.

For the filling capacity of the temperature control unit, refer to the data sheet. \rightarrow From page 94, section »Annex«. Or to the rating plate on the back of the temperature control unit. Please also consider: \rightarrow Page 27, section **»Ambient conditions**« and \rightarrow Page 29, section **»Installation condi**tions«.

Classifying the

e applica- tion field	Class of application field	Application field	Example of the installation location	Max. quantity of refrigerant		Max. permissible quantity above ground level (GL)
	A	General	Publicly accessible area in a public building		AND	1.5 kg
	В	Monitored	Laboratories	8 g/m ³	1	2.5 kg
	С	Access only for authorized persons	Production equipment	ambient air		10.0 kg
	Temperature	control units with m	ore than 1 kg refrigerant must not	t he installed held	w group	d level (GL)

Temperature control units with **more than 1 kg** refrigerant **must not be installed below ground level** (GL).

Temperature control units with up to 150 g natural refrigerant


- The temperature control unit has been constructed to the requirements of EU and EFTA countries.
- Use the table as guidance for classifying the application field. Respect the max. refrigerant quantity stated therein.

Temperature control units WITH pre-installed gas warning sensor and > 150 g natural refrigerant

- The temperature control unit has been constructed to the requirements of EU and EFTA countries.
- Use the table as guidance for classifying the application field. Respect the max. refrigerant quantity or the permissible highest quantity above ground level (GL) stated therein.
- Ventilation with optional supply and exhaust air connection: Use the temperature control unit's supply and exhaust air connection to connect it to the building's exhaust system. For the exact position please refer to the wiring diagram. \rightarrow From page 94, section **»Annex**«. First, remove the cover to the air inlet connection; an air filter mat is installed behind it. This air filter mat must be checked / replaced at regular intervals so that the air flowing into the temperature control unit is not reduced. → Page 81, section **»Function check and visual inspection**«. Connect the building's exhaust system with the temperature control unit's exhaust air port. The cover of the supply air port must not be removed if the exhaust system provided in a building is not used.

- The mounting plate for mounting a gas warning sensor is located inside the temperature control unit in the vicinity of the >Cable entry gas warning sensor< [100].
- For the position of the >Cable entry gas warning sensor< [100] please refer to the wiring diagram.
 → From page 94, section »Annex«.
- Additional information on the pre-installed gas detection sensor:
 - The built-in gas detection sensor enables a safety shutdown at 20% of the lower explosive limit via a power disconnect relay that is to be installed by the operator. The temperature control unit is thus switched off early and safely in case of fault.
 - A 24 V DC external power supply must be available for the pre-installed gas warning sensor. The alarm output of the gas warning sensor uses a 4 20 mA signal. Please refer to the data sheet of the gas warning sensor for further technical information. A separate processing unit is available as an accessory for the control of the power disconnect relay. The processing unit provides a potential-free switching contact and simultaneously provides the power supply and analysis of the gas warning sensor. Both variants require the operator to provide the necessary dimensioning and installation. Please refer to the data sheet of the gas detection system can be connected to the operator's alarm control unit. The operator is responsible for this and for the other measures.
 - The operator is responsible for the calibration of the gas detection sensor prior to initial operation and the observance of calibration and maintenance intervals according to the operating manual. We recommend to set calibration and maintenance intervals between 6 and 12 months if no information is provided. For increased safety requirements, shorter intervals can be specified. On request we will recommend a specialist company to carry out the calibration and maintenance.

Temperature control units WITHOUT pre-installed gas warning sensor and > 150 g natural refrigerant

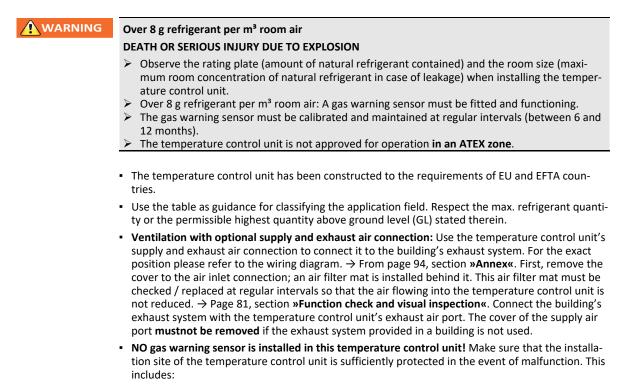

- Installation of a building's gas warning sensor (room monitoring).
- Permanent ventilation of the temperature control unit and/or the installation site.
- All-pole disconnection in the event of malfunction of the temperature control unit.

19

1.4.1.3 Temperature control units with fluorinated greenhouse gases/refrigerants

F gases regulation (EC) No. 517/2014 of April 16, 2014, on fluorinated greenhouse gases, and repealing Regulation (EC) No. 842/2006.

These regulations deal with all systems that contain fluorinated refrigerants. The substances dealt with in Directive (EC) No. 1005/2009 of the European Parliament and of the Council of 16 September 2009 that deplete the ozone layer are excluded (CFC/HCFC).

The directive regulates the reduction of the emission, utilization, recovery, and destruction of certain fluorinated greenhouse gases. It also regulates the identification and disposal of products and devices that contain these gases. Since July 4, 2007, responsible bodies must check their stationary refrigeration systems for leaks at regular intervals, and have any leaks eliminated immediately.

Directive (EC) No. 303/2008 contains stipulations on the training and certification of companies and personnel that are permitted to execute the specified activities.

Obligations of the responsible bodies:

- Directive (EC) No. 842/2006 already imposed a number of obligations upon responsible bodies regarding certain fluorinated greenhouse gases. The new Ordinance on Fluorinated Greenhouse Gases upholds these to a large extent. Some duties are added while others are designed differently by this new ordinance. Please refer to the text of this ordinance for a complete overview of the individual responsibilities of responsible bodies.
- General obligation to reduce emissions.
- Only certified companies may maintain, repair or decommission refrigeration systems. The responsible bodies must verify that these companies are certified.
- Regular leak tests of stationary refrigeration systems by certified personnel (such as Huber service engineers). The required test interval is based on the refrigerant filling capacity and the type of refrigerant, converted to CO₂2 equivalent.
- Responsibility of responsible bodies operating a plant to recover F-gases by certified personnel.
- Obligatory documentation requirement in the refrigeration system's operation manual, specifying type and volume of refrigerant used or recovered. The responsible body must keep the records for at least 5 years after their creation and present it to the responsible authority upon request.
- Temperature control units with natural refrigerants (NR) are exempt from this Directive.
- Please refer to the data sheet or rating plate of your temperature control unit for the quantity and type of refrigerant.
- Additional information about the definition of the inspection interval can be found on our website.

1.4.2 Requirements for operators

Work on the temperature control unit is reserved for appropriately qualified specialists, who have been assigned and trained by the responsible body to do so. Operators must be at least 18 years old. Under 18-year olds may operate the temperature control unit only under the supervision of a qualified specialist. The operator is responsible vis-a-vis third-parties in the work area.

1.4.3 Obligations of the operators

Carefully read the operation manual before operating the temperature control unit. Please observe the safety instructions. When operating the temperature control unit, wear appropriate personal protective equipment (e.g. safety goggles, protective gloves, non-slip shoes).

1.5 General information

1.5.1 Description of workstation

The workstation is located at the control panel in front of the temperature control unit. The workstation is determined by the customer's connected peripheries. Accordingly, it must be designed safe by the responsible body. The workstation design also depends on the applicable requirements of the German occupational health and safety regulations [BetrSichV] and the risk analysis for the workstation.



1.5.2 Safety devices to DIN 12876

The rating of your temperature control unit is stated on the data sheet in the appendix.

Rating of laboratory thermostats and	Classification	Temperature control medium	Technical requirements	Identification ^{d)}		
laboratory baths	I	Non-combustible ^{a)}	Overheat protection ^{c)}	NFL		
	Ш	Combustible ^{b)}	Adjustable overheat protection	FL		
	Ш	Combustible ^{b)}	Adjustable overtemperature protection and additional low-level protection	FL		
	^{b)} The temperature con ^{c)} The overheat protect	 ^{a)} Usually water; other fluids only if non-combustible even within the temperature range of an individual fault. ^{b)} The temperature control media must have a fire point of ≥ 65 °C. ^{c)} The overheat protection can, for instance, be realized using a suitable fill level sensor or a suitable temperature limiter. ^{d)} Optional at the choice of the manufacturer. 				
Overview of the temperature thresholds		Over-temperature tresho 5 K below combustion po Aaximum setpoint] oper threshold for temper ietpoint] an only be set in defined etween maximum/minir Ainimum setpoint] ower thershold for temper Ainimum working temper	ure of the temperature control unit. Id] int of the thermal fluid. erature control set by responsible body. range mum setpoint). erature control set by responsible body.			

The workings of the overtemperature and low-level protection depend the type of temperature control unit.

1.5.2.1 Temperature control unit (without heater) with classic float switch

Probably the most frequently found type today is the level monitor with a mechanical float switch. In the bath vessel, a floating body, which is guided in a device, floats on the surface of the thermal fluid. Depending on the level of the thermal fluid, the float device signals the electronics a **state of good** (in case of sufficient filling) or a **state of bad** (in case of insufficient filling). You should check the functionality of the float device from time to time. To do so and when in stand-by mode, press the float body into the bath with a tool (e.g. a screwdriver). The electronics must report an alarm.

1.5.2.2 Temperature control units (with heater) with electronic low level sensor

Combined electronic low level and overtemperature protection

Some temperature control units (device dependent) feature an electronic overtemperature and low-level protection. Instead of a mechanical float device, temperature sensors are employed on the surface of the heating coils. These sensors monitor overtemperature at this potential ignition source, thus ensuring that the controller regulates the heating coil temperatures below the critical temperature (fire point of the thermal fluid). This does not apply to chillers where the level is capacitively determined). An appropriate message is output on the **>Touchscreen< [88]**. Please note the figure **"Pilot ONE"**. \rightarrow Page 6.

There is no longer a need for mechanical tools to set the trigger values of the overtemperature protection. It is replaced by a software engineering tool. The threshold value for the overtemperature protection can be set only if a code, randomly generated by "Pilot ONE", is entered correctly. As with the mechanical tool, accidental settings are thus prevented.



1.5.3 Further protective devices

INFORMATION

Emergency strategy – interrupt the power grid connection! To determine the type of switch or switch combination your temperature control unit is equipped with, please refer to the wiring diagram. \rightarrow From page 94, section **»Annex**«.

Temperature control units with >Mains isolator< [36] (red/yellow or gray): Turn the >Mains isolator< [36] to the "0" position.

Temperature control units with >Mains isolator< [36] (red/yellow) and additional >Appliance switch< [37] (gray): Turn the >Mains isolator< [36] to the "0" position. Then turn the >Appliance switch< [37] to the "0" position!

Temperature control units with >Mains isolator< [36] (gray) and >Emergency stop switch< [70] (red/yellow): Press the >Emergency stop switch< [70]. Then turn the >Main switch< [36] to the "0" position!

Temperature control units with >Mains switch< [37]: Power supply via socket: Disconnect the temperature control unit from the power supply. Then turn the **>Mains isolator< [37]** to the "0" position! Power supply via hard wiring: Disconnect the power grid supply by means of the building's circuit breaker. Then turn the **>Mains isolator< [37]** to the "0" position!

Temperature control units without a switch or inside a protective housing: Connection via socket: Disconnect the temperature control unit from the power supply. Connection via hard wiring: Disconnect the power grid supply by means of the building's circuit breaker!

1.5.3.1 Power interruption

Following a power outage (or when switching on the temperature control unit), this function can be used to determine how the temperature control unit is supposed to respond. This response can be determined via "Pilot ONE".

OFF/Standby (Default setting)

After turning the temperature control unit on, thermoregulation is started only after manual input.

ON / Temperature control active

After turning on the temperature control unit, thermoregulation is always started. An INFO appears for a few seconds. This makes it possible to suppress the automatic start.

Power Failure Auto

If temperature control is active during a power outage, the process will automatically continue after the power outage.

1.5.3.2 Alarm functions


An alarm is a system state that signals unfavorable process conditions. The temperature control unit can be programmed so that the plant operator is warned when defined limit values are exceeded. The response of the temperature control unit to an alarm can be determined. Possible responses are: Switch off temperature control or control temperature to a safe setpoint (2nd setpoint).

1.5.3.3 Warning messages


Warning messages contain a message about the irregularity of the temperature control unit. These messages have no further consequences. The plant operator evaluates the relevance of the message and takes action where necessary.

Example: Air and water cooling

Chapter 1

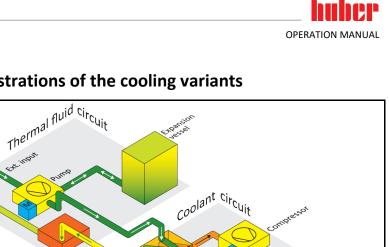



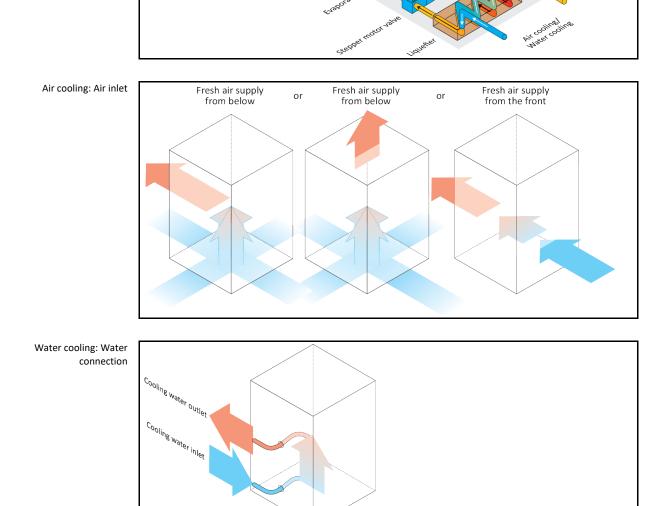

cc* 23



1.6.1 Consequence of inadequate energy dissipation

Room air/cooling water

Consequences of, for instance, contamination of the liquefier fins, inadequate clearance between temperature control unit to wall/bath wall, room air/cooling water too warm, cooling water differential pressure too low, suction strainer contamination: The refrigerant in the coolant circuit can no longer fully discharge the admitted energy to the room air/cooling water. Thus there is not sufficient liquefied refrigerant available, the condensation temperature and the energy consumption to rise.

Coolant circuit

Consequences of inadequate refrigerant quantity/rising condensation temperature: Not all the cooling capacity from the coolant circuit is available at the evaporator. This means reduced energy transmission from the thermal fluid circuit.

Thermal fluid circuit

Consequence of inadequate energy dissipation from the thermal fluid: The thermal fluid can only dissipate the energy from your application to a limited extent.

Application

Consequences of inadequate energy dissipation from the application: The energy created (exothermic) in the application can no longer be fully dissipated.

Temperature control unit

An electronically-controlled expansion valve is used in the temperature control unit to optimize the power adjustment. The expansion valve always provisions the maximum possible cooling capacity within the permissible ambient temperature range. The temperature control unit switches off when the upper range is reached (maximum permissible ambient temperature).



2 Commissioning

2.1 In-plant transport

WARNING	Temperature control unit is not transported / moved according to the specifications in this oper- ation manual		
	DEATH OR SERIOUS INJURY DUE TO CRUSHING		
	 Always transport / move the temperature control unit according to the specifications in this operation manual. Wear personal protective equipment during transport. Always work with the specified number of persons when moving the temperature control unit on casters (if any). If the temperature control unit is equipped with casters and parking brakes: 2 parking brakes are always freely accessible when moving the temperature control unit. Activate the 2 parking brakes in an emergency! If only one parking brake is activated on the casters in an emergency: The temperature control unit is not stopped but rotates around the axis of the caster with the activated parking brake! 		
ΝΟΤΕ	Temperature control unit transported in a horizontal position		
NOTE	DAMAGE TO THE COMPRESSOR		
	 Only transport the temperature control unit in an upright position. 		
ΝΟΤΕ	Filled temperature control unit is transported		
	MATERIAL DAMAGE DUE TO OVERFLOWING THERMAL FLUID		
	 Only transport an emptied temperature control unit. 		
	 If available, use the lugs on the top side of the temperature control unit for transportation. Use an industrial truck for transport. The casters (if present) on the temperature control unit are not suitable for transport. The casters are symmetrically loaded with 25% of the total mass of the temperature control unit. Remove the packing material (e.g. the palette) only at the place of installation. Protect the temperature control unit from transport damage. Do not transport the temperature control unit alone and without aids. Check the load bearing capacity of the transportation route and the place of installation. The parking brakes at the casters (if any) must be activated and/or the leveling feet (if any) must be unscrewed/activated before the temperature control unit is put into operation. → Page 33, section »Unscrewing/activating the leveling feet (if any)«. 		
2.1.1	Lifting and transporting the temperature control unit		
2.1.1.1	Temperature control unit with lifting eyes		

 NOTE
 The temperature control unit is raised at the lifting eyes without load handling attachments

 DAMAGE TO THE TEMPERATURE CONTROL UNIT
 Always use load handling attachments when lifting and transporting the temperature control unit.


 > The lifting eyes are only designed for a load without inclination (0°).
 The load handling attachment used must be adequately dimensioned. Take the dimensions and weight of the temperature control unit into account.

 Example: lifting eyes (round, angular, and recessed (left to right))
 Image: Description of the temperature control unit into account.

- Do not lift and transport the temperature control unit at the lifting eyes alone and without aids.
- Lift and transport the temperature control unit at the lifting eyes only with a crane or an industrial truck.
- The crane or industrial truck must have a lifting force equal to or greater than the weight of the temperature control unit. See the data sheet for the weight of the temperature control unit.
 → From page 94, section »Annex«.
- If the leveling feet have been removed for shipping: Only lower the temperature control unit when all leveling feet have been installed. → Page 26, section »Mounting/removing leveling feet«.

2.1.1.2 Temperature control unit without lifting eyes


Example: Supporting points for forklift arms for free-standing models from a certain overall size. For the exact position please refer to the wiring diagram in the annex.

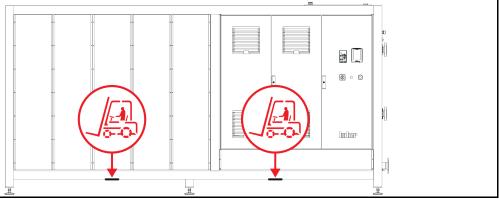

- Do not lift and transport the temperature control unit alone and without aids.
- · Lift and transport the temperature control unit only with an industrial truck.
- The industrial truck must have a lifting force equal to or greater than the weight of the temperature control unit. See the data sheet for the weight of the temperature control unit. → From page 94, section »Annex«.
- If the leveling feet have been removed for shipping: Only lower the temperature control unit when all leveling feet have been installed. → Page 26, section »Mounting/removing leveling feet«.

2.1.2 Mounting/removing leveling feet

Only valid if the leveling feet have been removed for shipping.

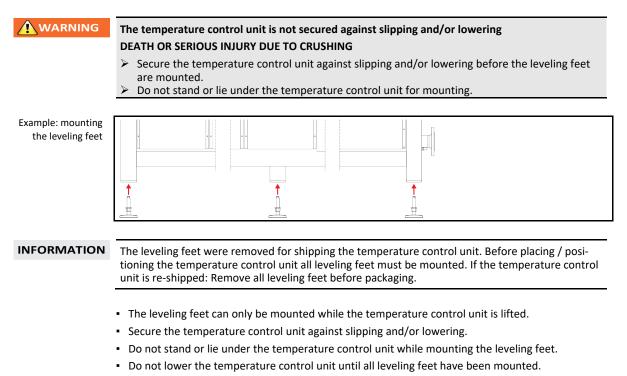



2.1.3 Positioning the temperature control unit

2.1.3.1 Temperature control unit with casters

- Do not use the casters for the transportation to the place of installation. → Page 25, section »Lifting and transporting the temperature control unit«.
- Use the casters only for positioning at the place of installation.
- Only ever move the temperature control unit on the casters if the surface is level, without gradient, non-slip and stable.
- Do not move the temperature control unit alone.
- At least 2 persons are required to move the temperature control unit on casters. At least 5 persons are required to move the temperature control unit on the casters if the total weight of the temperature control unit is over 1.5 tons.
- The parking brakes must be activated at the casters and/or the feet (if present) must be unscrewed/activated before the temperature control unit is put into operation. → Page 33, section »Unscrewing/activating the leveling feet (if any)«.

2.1.3.2 Temperature control unit without casters

- An industrial truck must be used for positioning the temperature control unit.
- Do not move the temperature control unit alone.
- At least 2 persons are required to move the temperature control unit.
- The industrial truck must have a lifting force equal to or greater than the weight of the temperature control unit. See the data sheet for the weight of the temperature control unit. → From page 94, section »Annex«.
- The leveling feet (if present) must be unscrewed/activated before the temperature control unit is
 put into operation.
 → Page 33, section »Unscrewing/activating the leveling feet (if any)«.

2.2 Unpacking

WARNING

Starting up a damaged temperature control unit DANGER TO LIFE FROM ELECTRIC SHOCK

- > Do not operate a damaged temperature control unit.
- > Please contact Customer Support. \rightarrow Page 92, section **»Contact data**«.

PROCEDURE

- Check for damage to the packaging. Damage can indicate material damage to the temperature control unit.
- > Check for any transport damage when unpacking the temperature control unit.
- Always contact your forwarding agent regarding the settlement of claims.
- ➢ Observe the proper disposal of packaging material. → Page 17, section »Proper disposal of resources and consumables«.


2.3 Ambient conditions

 ▲ CAUTION
 Unsuitable ambient conditions / unsuitable installation

 SERIOUS INJURY DUE TO CRUSHING
 > Comply with all requirements! → Page 27, section »Ambient conditions« and → Page 29, section »Installation conditions«.

 INFORMATION
 Make sure there is adequate fresh air available at the site for the circulation pump and the compressors. The warm exhaust air must be able to escape upwards unhindered.

 Free-standing models
 For the connection data, see the data sheet. → From page 94, section »Annex«.

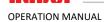

Use of the temperature control unit is permitted only under normal ambient conditions in accordance with the currently valid DIN EN 61010-1.

- Use only indoors. The illuminance must be at least 300 lx.
- Installation altitude up to 2,000 meters above sea level.
- Maintain wall and ceiling clearance for adequate air exchange (dissipation of waste heat, supply
 of fresh air for the temperature control unit and work area). Ensure adequate floor clearance for
 air-cooled temperature control units. Do not operate this temperature control unit from within
 the box or with an inadequately dimensioned bath. This inhibits the air exchange.
- Ambient temperature values are provided on the technical data sheet; to ensure trouble-free operation, compliance with the ambient conditions is mandatory.
- Relative humidity max 80% to 32 °C and 40 °C decreasing linearly to 50%.
- Short distance to supply connections.
- The temperature control unit must not be installed so as to hinder or even prevent access to the disconnecting device (to the power supply).
- For the magnitude of the mains voltage fluctuations, refer to the datasheet. → From page 94, section »Annex«.
- Transient surges, as would normally occur in the power supply system.
- Installation Class 3
- Applicable degree of soiling: 2.
- Surge category II.

Please note: \rightarrow Page 23, section **»Exemplary illustrations of the cooling variants«**.

Wall clearances		Distance in cm		
		Air cooling	Water cooling	
	Side			
	[A1] Top	Air outlet on top of unit: free standing	-	
	[A2] Top	can be located under a bench	can be located under a bench	
	[B] Left	min. 20	min. 10	
	[C] Right	min. 20	min. 10	
	[D] Front	min. 20	min. 10	
	[E] Rear	min. 20	min. 20	



	Distance in cm (for	operation in a tub)	
	Air cooling	Water cooling	
Side	[A1] $[B]$ $[C]$ $[D]$ $[E]$ $[E]$		
[A1] Top	Air outlet on top of unit: free standing	_	
[A2] Top	can be located under a bench	can be located under a bench	
[B] Left	min. 20	min. 20	
[C] Right	min. 20	min. 20	
[D] Front	min. 20	min. 20	
[E] Rear	min. 20	min. 20	

2.3.1 EMC-specific notes

INFORMATION Connecting cables in general

Prerequisites for a failure-free operation of the temperature control units incl. their connections with external applications: Installation and wiring must be carried out professionally. Related topics: "Electrical safety" and "EMC-compliant wiring".

Cable lengths

For flexible/fixed cable routing of more than 3 meters, the following must amongst other things be observed:

 Equipotential bonding, grounding (see also technical data sheet "Electromagnetic compatibility EMC")

- Compliance with "external" and/or "internal" lightning/overvoltage protection.

- Design protection measures, professional cable selection (UV resistance, steel pipe protection, etc.)

Attention:

The operating company is responsible for compliance with national/international directives and laws. This also includes the testing of the installation/wiring required by law or standards.

This device is suitable for operation in **"industrial electromagnetic environments"**. It meets the **"immunity requirements"** of the currently applicable **EN61326-1**, which are required for this environment.

It also meets the **"interference emission requirements"** for this environment. It is a **Group 1** and **Class A** device according to the currently applicable **EN55011**.

Group 1 specifies that high frequency (HF) is only used for the function of the device. **Class A** defines the interference emission limits to be observed.

2.4 Installation conditions

WARNING

Temperature control unit is connected to the power supply line DEATH FROM ELECTRICAL SHOCK BY DAMAGE TO THE POWER CABLE.

Do not put temperature control unit on power cable.



Operating the temperature control unit fitted with casters without brakes activated

CRUSHING OF LIMBS

Activate brakes on the casters.

- Allow the temperature control unit to acclimate for about 2 hours when changing from a cold to a warm environment (or vice versa). Do not turn on the temperature control unit before!
- Install upright, stable and tilt-resistant.
- Use a non-combustible, sealed foundation.
- Keep the environment clean: Prevent slip and trip hazards.
- Wheels, if installed, must be locked after installation!
- Spilled/leaked thermal fluid must be removed immediately. Observe the proper disposal of thermal fluid and aids. → Page 17, section »Proper disposal of resources and consumables«.
- Observe the floor load bearing capacity for large units.
- Observe the ambient conditions.

2.5 Recommended temperature control and cooling water hoses

	Use of unsuitable/defective hoses and/or hose connections		
	INJURIES		
	> Thermal fluid		
	Use appropriate hoses and/or hose connections.		
	Check periodically for leaks and the quality of the hose and hose connections and take suitable		
	 measures (replace) as required. Isolate and protect temperature control hoses against contact/mechanical load. 		
	 Cooling water 		
	 Reinforced hoses must be used to satisfy tougher safety requirements. 		
	Shut off the cooling water supply to the temperature control unit even for shorter downtimes		
	(e.g. overnight).		
	Hot or cold thermal fluid and surfaces		
	BURNS TO LIMBS		
	Avoid direct contact with the thermal fluids or the surfaces.		
	Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).		
	guggies, salety luotweal).		
	Uncontrolled formation of ice at the connections and hoses of the thermal fluid circuit		
<u></u>	SLIP AND TRIP HAZARD		
	 If the temperature is controlled in the minus range, ice forms at the hoses and connections of 		
	the thermal fluid circuit. This occurs by condensing and freezing of atmospheric humidity.		
	> Check the strength of the ice formation. If too much ice is formed, this increases the risk of the		
	temperature control unit tipping over. Secure the temperature control unit against tipping if		
	this is the case.Check the ground below the ice formation for condensation water. Collect the condensation		
	water with a suitable container or thoroughly remove it at regular intervals. You thus prevent		
	the danger of slipping caused by condensation.		
	To connect applications, use only temperature control hoses that are compatible with the thermal		
	fluid used. When selecting temperature control hoses, also pay attention to the temperature range		
	in which the hoses are to be used.		
	 We recommend you use only temperature-insulated temperature control hoses with your tem- perature control unit. The responsible body is responsible for the insulation of connection valves. 		
	 We exclusively recommend reinforced hoses for connecting to the cooling water supply. Cooling water and insulated temperature control hoses can be found in the Huber catalogue under Acces- 		

sories.



. .

2.6 Wrench sizes and torques

.

Note the wrench sizes that result for the pump connection on the temperature control unit. The following table lists the pump connections and the resulting wrench sizes, and torque values. A leak test must always be performed, and the connections tightened if necessary. The values of the maximum torque (see table) must **not** be exceeded.

.

. .


Overview wrench sizes and

Overview sizes and	Connection	Sleeve nut wrench size	Connector wrench size	Recommended torques in Nm	Maximum torques in Nm
torques	M16x1	19	17	20	24
	M24x1.5	27	27	47	56
	M30x1.5	36	32	79	93
		36	36	79	93
	M38x1.5	46	46	130	153
	G-thread (flat- sealing)	Adapt the torque to the material of the flat seal used. First hand-tighten the temperature control hose. When using adapters, do not overtighten the G-thread on the pump connection when connecting a temperature control hose. When connecting a temperature control hose to the adapter piece, secure the G thread against overwinding.			

2.7 Temperature control units with water cooling

WARNING	Open electrical wires below the temperature control unit if the temperature falls below the dew point.		
	DEATH FROM ELECTRICAL SHOCK BY WATER ENTRY INTO THE ELECTRIC LINES.		
	A temperature below the dew point may result in condensation in the temperature control unit and at the cooling water connections. The condensation is caused by high humidity at the cool- ing water-bearing components. The condensation exists the temperature control unit at the bottom. Electrical lines directly below the temperature control unit must be protocted expired lines.		
	Electrical lines directly below the temperature control unit must be protected against liquid ingress.		
	Steam or hot water escapes at the >cooling coil< [29]		
	SEVERE BURNS		
	Do not use the >cooling coil< [29] at bath temperatures > 95 °C.		
	Provide a greater flow rate through the >cooling coil< [29] at bath temperatures > 60 °C.		
	Use of unsuitable/defective hoses and/or hose connections		
	INJURIES		
	➢ Thermal fluid		
	 Use appropriate hoses and/or hose connections. 		
	Check periodically for leaks and the quality of the hose and hose connections and take suitable		
	measures (replace) as required.		
	Isolate and protect temperature control hoses against contact/mechanical load.		
	Cooling water Deinforced bases must be used to satisfy taugher safety requirements		
	 Reinforced hoses must be used to satisfy tougher safety requirements. Shut off the cooling water supply to the temperature control unit even for shorter downtimes 		
	(e.g. overnight).		
	(0.8.000111810).		
NOTE	No protection against corrosion		
	DAMAGE TO THE TEMPERATURE CONTROL UNIT		
	The addition of anti-corrosion agents is mandatory if salts (chlorides, bromide) have been		
	added to the water circuit.		
	Ensure that the materials used in the cooling water circuit are resistant with respect to the cooling		
	water. For information on materials used see the data sheet. \rightarrow From page 94, section »Annex« .		
	Take suitable measures to maintain the warranty conditions.		
	For information about water quality, see www.huber-online.com.		



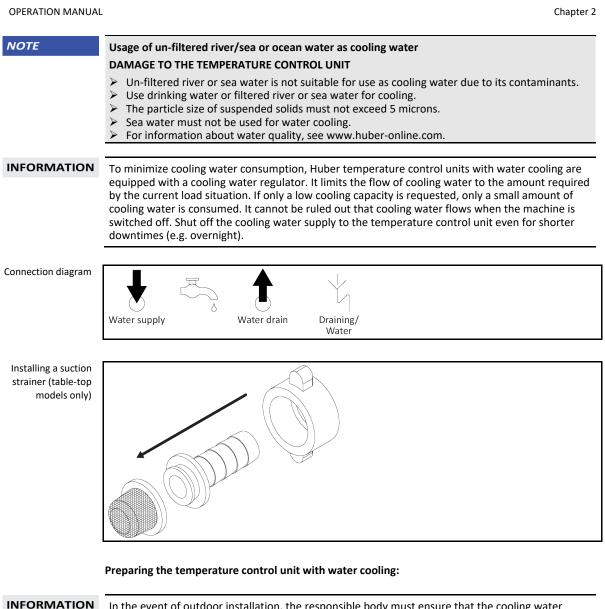

ON In the event of outdoor installation, the responsible body must ensure that the cooling water supply and return lines are laid frost-protected. The cooling water temperature must not fall below 3 °C. At ambient temperatures below 3 °C, the cooling water supply must be heated.

The minimum pressure differential in the cooling water circuit and the recommended cooling water inlet temperature can be found on the data sheet. \rightarrow From page 94, section **»Annex**«.

Observe the wiring diagram. \rightarrow From page 94, section **»Annex**«.

PROCEDURE

- Close (if fitted) the >Cooling water drain< [15].</p>
- > Connect the >Cooling water outlet< [14] to the water return flow. A seal must be used.
- Insert the suction strainer (dirt trap) into the >Cooling water return< [13].</p>
- Connect the >Cooling water inlet< [13] to the water supply.</p>

NOTE

Leaking cooling water connections DAMAGE BY ROOM FLOODING

- Slowly open the building-side shut-off valves of the cooling water supply and return line.
- If water leaks from the cooling water connections: shut off the cooling water supply and return line immediately.
- Provide leakproof cooling water connections.




- Open the shut-off valves in the water line on the temperature control unit (if present) and on the building side.
- Check the connections for leaks.

2.8 Temperature control unit with batch

2.8.1 Operation as bath thermostat

Figure bath thermo-

stat

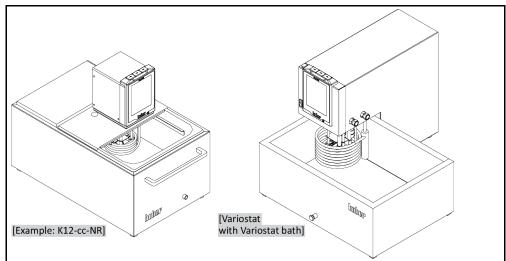

Note the volume displacement caused by a sample (e.g. Erlenmeyer flask). Place your sample into the empty bath. Only then fill in a sufficient amount of temperature control medium. Also note that the level of the temperature control medium drops when you remove the sample. This may cause a safety shutdown (low level protection) during an enabled thermoregulation. Therefore, switch off the temperature control unit beforehand.

The entire area/volume of the "bath vessel Variostat" is not available (see picture above) in the combination Variostat with "bath vessel Variostat". For bath thermostats: Please note the size of the bath opening. \rightarrow From page 94, section **»Annex«**.

2.9 Preparations for operation

2.9.1 Unscrewing/activating the leveling feet (if any)

The leveling feet are not unscrewed/activated before switching on the temperature control unit DEATH OR SERIOUS INJURY DUE TO CRUSHING

- The parking brakes must be activated at the casters (if any) and/or the leveling feet must be unscrewed/activated before the temperature control unit is put into operation.
- The temperature control unit may move if the parking brakes of the casters (if any) are not activated and/or the leveling feet are not unscrewed/activated.


Always unscrew/activate the leveling feet before switching on the temperature control unit. Uneven floors can be compensated by adjusting these leveling feet.

PROCEDURE

- > Verify that the parking brakes of the casters (if any) have been activated.
- Unscrew the leveling feet.
- Compensate uneven floors by adjusting these leveling feet, if necessary. Use a spirit level to horizontally align the temperature control unit.
- Tighten the lock screws on the leveling feet after aligning the temperature control unit. This prevents the leveling feet from changing their height during operation.



Chapter 2

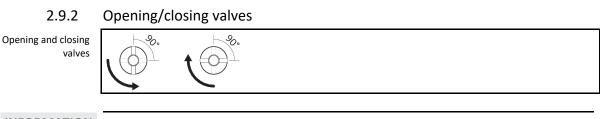

INFORMATION

Opening valves: Open valves by turning them counterclockwise (turn 90° left as far as it will go). **Close valves:**

close valves by turning them clockwise (turn 90° right as far as it will go).

PROCEDURE


- > Check that all valves are closed.
- Close all valves by turning them clockwise (turn 90° right as far as it will go).

2.9.3 Installing collecting container

PROCEDURE

- Install a suitable hose at the >overflow< [12] on the temperature control unit (if present). The hose must be compatible with the thermal fluid and the temperature.</p>
- > Put the other end of the hose in a suitable collecting container.

2.9.4 Cooling/heating baths: Connect the energizing lead

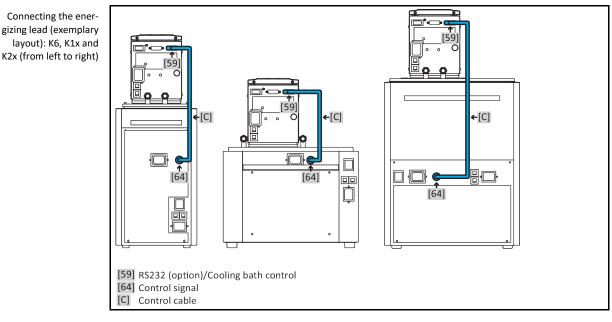

INFORMATION

The energizing lead must be installed only if the cooling bath is used **in combination** with an immersion circulator thermostat.

PROCEDURE

Connect the >RS232 (option)/activation cooling bath< [59] port on the hanger thermostat with the >activation connector< [64] on the cooling bath. The cable required is included.</p>

NOTE



2.9.5 Connecting the functional earth

PROCEDURE

➢ If required, connect the temperature control unit's >Functional ground terminal< [87] to the building's grounding point. Use a ground strap for this purpose. For the exact position and thread size please refer to the wiring diagram. → From page 94, section »Annex«.</p>

2.10 Connecting externally closed application

Observe the wiring diagram. \rightarrow From page 94, section **»Annex**«.


2.10.1 Connecting an externally closed application

Overpressure in the application (e.g. > 0.5 bar (g) with glass apparatus) MATERIAL DAMAGE TO THE APPLICATION

- Provide an overpressure protective device to prevent damage to the application.
- Do not install valves/quick-release couplings in the feed/discharge lines from the temperature control unit to the application and from the application to the temperature control unit.
- If valves/quick-release couplings are required:
- Install burst disks on the application itself (at the feed and discharge lines).
- Install a bypass upstream of the valves/quick-release couplings for the application.
 Matching accessories (e.g. bypasses to reduce pressure) can be found in the Huber catalog.

Example: Connecting an externally closed

application

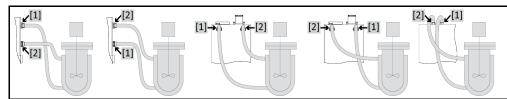

To enable your application to be operated correctly and eliminate air bubbles from the system, you must ensure that the **>Circulation flow< [1]** connection from the temperature control unit is attached to the lower connection point of the application and the **>Circulation return< [2]** into the temperature control unit is attached to the higher connection point of the application.

PROCEDURE

- Remove the screw plugs from the >Circulation flow< 1 and >Circulation return< 2 connections.</p>
- Then connect your application to the temperature control unit using suitable thermal fluid hoses. Observe the table with the wrench sizes. \rightarrow Page 31, section **»Wrench sizes and torgues**«.
- Observe the table with the wrench sizes. → Page 31, section **»wrench sizes and to**Check the connections for leaks.

2.11 Connecting to the power supply

INFORMATION Based on local circumstances, it may be that you need to use an alternative power cable instead of the supplied original power cable. Do not use a power cable that is longer than **3 m** to be able to disconnect the temperature control unit at any time from the mains. Have the mains cable only replaced by a qualified electrician.

2.11.1 Connection using socket with protective earth (PE)

Connecting to a power socket without protective earth (PE) MORTAL DANGER FROM ELECTRIC SHOCK

> Always connect the temperature control unit to safety sockets (PE).



OPERATION MANUA	L Chapter 2
DANGER	 Damaged power cable/power cable connection MORTAL DANGER FROM ELECTRIC SHOCK Do not start up the temperature control unit. Isolate the temperature control unit from the power supply. Have the power supply cable/power supply connection replaced and inspected by an electrician. Do not use a power cable that is longer than 3 m.
ΝΟΤΕ	 Incorrect power supply connection DAMAGE TO THE TEMPERATURE CONTROL UNIT Your building's existing power supply voltage and frequency must match the data provided on the rating plate of the temperature control unit.
INFORMATION	In case of uncertainties about an existing protective earth (PE), have the connection inspected by an electrician.
2.11.2	Connection via hard wiring
DANGER	Connection/adjustment to the power supply not carried out by an electrician MORTAL DANGER FROM ELECTRIC SHOCK ➤ Have the connection/adjustment to the power supply carried out by an electrician.
DANGER	 Damaged power cable/power cable connection MORTAL DANGER FROM ELECTRIC SHOCK Do not start up the temperature control unit. Isolate the temperature control unit from the power supply. Have the power supply cable/power supply connection replaced and inspected by an electrician. Do not use a power cable that is longer than 3 m.
ΝΟΤΕ	 Incorrect power supply connection DAMAGE TO THE TEMPERATURE CONTROL UNIT Your building's existing power supply voltage and frequency must match the data provided on the rating plate of the temperature control unit.
2.11.3	Converting the power supply connection
INFORMATION	To restore the factory settings, the connection values must be re-entered .
	To find out whether your temperature control unit is convertible from 230V 1 $^{\sim}$ 50Hz to 400V 3 $^{\sim}$ N 60Hz, refer to the datasheet. \rightarrow From page 94, section »Annex «.
	 The maximum current consumption at 230 V operation is just below 16 A. However, there are 230 V power supplies with a smaller fuse protection. During initial startup you will be prompted to adjust the power consumption of the temperature control unit to your power supply. The following can be selected: 16 A (100 % heating power with switched on compressor) 13 A (reduced heating power with switched on compressor) 10 A (reduced heating power with switched on compressor) A disconnected compressor does not effect the heating power: 16 A (100 % heating power with switched off compressor) 13 A (100 % heating power with switched off compressor) 13 A (100 % heating power with switched off compressor) 13 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power with switched off compressor) 10 A (100 % heating power heating power with switched off compressor) 10 A (100 % he

Installing a mains connection (exemplary layout) - pure cooling bath (K6, K1x and K2x (from left to right), valid for 100 V, 115 V

and 230 V versions)

Chapter 2

2.11.4 Cooling/heating baths: Connecting to the power supply

Image: Second system Image: Second

2.11.4.1 Utilization as a pure cooling bath (without hanger thermostat)

PROCEDURE

Connect the cooling bath to the power supplies (valid for 100 V, 115 V and 230 V versions)


Connect the >power supply< [35] on the cooling bath with the building's power supply. The cable required is included.</p>

2.11.4.2 Utilization as a cooling/heating baths (with hanger thermostat)

There are **two** possibilities for the connection with the building's power supply. **1. Connection with two separate power supplies (valid for 100 V, 115 V and 230 V versions)** The cold bath and the hanger thermostat are **each separately** connected to the building's power supply.

2. Connection with only one power supply (only valid for 230 V version)

Only the cooling bath is connected to the building's power supply, while the hanger thermostat is connected to the cooling bath for a supply of power.

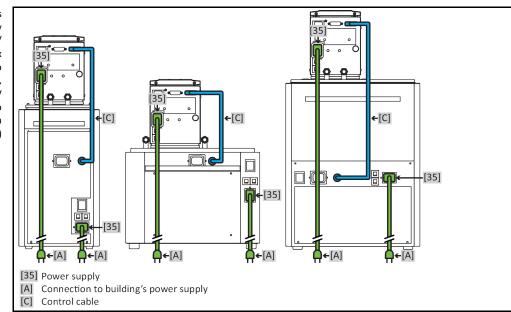

Installing a mains connection (exemplary layout), cooling bath / heating bath (K6, K1x and K2x (from left to right), valid for 100 V, 115 V and 230 V versions with two power supplies from the building)

INFORMATION




PROCEDURE

Connect the cooling bath and the hanger thermostat with TWO separate power supplies (valid for 100 V, 115 V and 230 V versions)

- Connect the >power supply< [35] on the hanger thermostat with the building's power supply. The cable required is included.
- Connect the >power supply< [35] on the cooling bath with the building's power supply. The cable required is included.</p>

Installing a mains connection (exemplary layout), cooling bath / heating bath (K6, K1x and K2x (from left to right), valid for 230 V versions with a single power supply from the building)

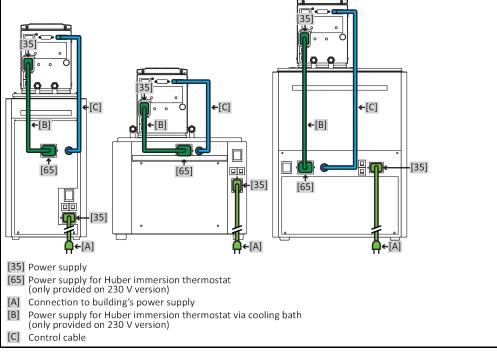

PROCEDURE

Connection of the cooling bath and the hanger thermostat with ONE power supply (only valid for 230 V version)

- Connect the >power supply< [35] on the hanger thermostat with the >power supply< [65] on the cooling bath. The cable required is included.</p>
- Connect the >power supply< 35 on the cooling bath with the building's power supply. The cable required is included.</p>



3 Function description

3.1 Function description of the temperature control unit

3.1.1 General functions

This temperature control unit is designed to be used with the **internal bath** as well as **externally closed applications**. \rightarrow Page 35, section **»Connecting an externally closed application«**. This **controller generation** is characterized by the fact that there is only a **single hardware**. For a license fee you will receive an **activation key**, which allows you to quickly expand the functionality of e.g. a simple temperature control unit (basic version) to its **maximum functionality** (professional version) at any time.

Thanks to the **low internal volume** and the **efficient refrigeration/heat technology** (in case of temperature control units with compressor cooling, you will achieve relatively **short cooling/heating rates**.

The cooling baths are designed to be used as both **pure cooling baths** as well as **in combination** with an **immersion circulator thermostat** (cooling/heating baths). In combination with an immersion circulator thermostat, the cooling baths may be used across the entire specified temperature range and can cool at the maximum operating temperature in continuous operation too.

3.1.2 Other functions

The **speed-controlled pump**, integrated into many table-top models and some free-standing models, enables the **speed** and thus the **pressure** of the thermal fluid to be controlled and hence optimally matched to the specified application.

The **self-optimizing cascade control** delivers **optimal controller results** for regulating levels up and down, both with **setpoint changes** and also with **exothermic reactions**. Temperature control can be either a-periodic or with a slight overshoot (faster).

Your temperature control unit can be easily integrated into many laboratory automation systems. This is provided for by the existing **standard interfaces Ethernet**, **USB device and USB host at the** "**Pilot ONE**". The optional Com.G@te allows you to expand your temperature control unit by digital interfaces (RS232 and RS485), an analogue current loop interface (0/4 - 20 mA or 0 - 10 V) and various digital control options (In/Out).

The **removable control panel ("Pilot ONE")** can also be used as a **remote control**. Please contact your dealer or Huber Sales Department if you need an extension cable. \rightarrow Page 92, section **»Contact data«**.

A connection jack for Pt100 process controller sensor enables you to accomplish external temperature control tasks with ease (requires E-grade Exclusiv or Professional).

The temperature control unit can be retrofitted with an **integrated temperature ramp function** and an **internal temperature programmer** via an "E-grade" upgrade. The "E-grade" upgrade **"Exclusiv"** enables you to set and call 3 temperature control programs with 5 program steps each by using the integrated programmer. The "E-grade" upgrade **"Professional"** enables you to set and call 10 temperature control programs with a total of 100 program steps by using the integrated programmer.

The temperature control unit has **overtemperature protection to DIN EN 61010-2-010, independent** of the control circuit itself.



Information on the thermal fluids 3.2

	 INJURIES Risk of injury to the eyes, sl The safety data sheet for the content must be respected. 	e thermal fluid to be used must be read prior to using it and its				
	 > Observe the local regulations/work instructions. > Wear your personal protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear). > Danger of slipping because floor and work area are contaminated. Clean the workplace; observe the proper disposal of thermal fluid and aids. → Page 17, section »Proper disposal of resources and consumables«. 					
ΝΟΤΕ	Non-compliance with the com unit	patibility between the thermal fluid and your temperature control				
	MATERIAL DAMAGE					
	 Observe the classification of your temperature control unit according to DIN 12876. Ensure the following materials are resistant with respect to the thermal fluid: Stainless steel 1.4301/ 1.4401 (V2A), copper, nickel, FKM, red bronze/brass, silver solder and plastic. The maximum viscosity of the thermal fluid must not exceed 50 mm²/s at the lowest working temperature! The maximum density of the thermal fluid may not exceed 1 kg/dm³! 					
ΝΟΤΕ	etc.) in a thermofluid circuit The thermal fluid circuit mu	of thermofluid (such as mineral oil, silicone oil, synthetic oil, water,				
Thermal fluid: Water	Designation	Specification				
	Calcium carbonate per liter	\leq 1.5 mmol/l; corresponds to a water hardness of: \leq 8.4 °dH (soft)				
	pH value	between 6.0 and 8.5				
	Ultrapure water, distillates	Add 0.1 g of sodium carbonate (Na $_2$ CO $_3$) per liter				
	Non-approved water	Distilled, deionized, demineralized, chloric, ferruginous, ammoniacal, con- taminated or untreated river water or sea water				
	Volume circulated (at least)	3 l/min. (not valid for cooling baths)				
	Thermal fluid: Water without ethylene glycol					
	Use	≥ +5 °C				
	Thermal fluid: Water-ethylene glycol mixture					
	Use	< +5 °C				
	Thermal fluid composition	The mixture's temperature must be 10 K below the permissible min. tem- perature. For the permissible temperature range, refer to the datasheet. → From page 94, section »Annex «.				

INFORMATION For thermal fluids we recommend the media listed in the Huber catalog. The name of a thermal fluid is derived from its working temperature range and its viscosity at 25 °C.



3.2.1 Only valid for CC[®]-520w, CC[®]-525w, CC[®]-906w and cooling baths used in continuous operation

Thermal fluid: Water

Designation	Specification			
Thermal fluid: Water without ethyl	ene glycol			
Use	do not use			
Thermal fluid: Water-ethylene glycol mixture				
Use	do not use			

3.3 To be noted when planning the test

INFORMATION Observe the intended operation. \rightarrow Page 15, section **»Proper operation**«. The focus is on your application. Bear in mind that system performance is influenced by heat transfer, temperature, thermal fluid viscosity, volume flow, and flow speed. • Make sure the electrical connection is adequately dimensioned. • The installation location of the temperature control unit should be selected so as to ensure adequate fresh air, even with water-cooled chillers. The maximum flow pressure of a temperature control unit must be taken into account in case of pressure-sensitive applications, such as glass reactors. Cross-section reduction or shut-off in the thermal fluid circuit must be avoided. Take appropriate measures to limit the pressure in the system. Observe the data sheet of the temperature control unit and the glass apparatus. \rightarrow From page 94, section **»Annex**«. Check whether it is necessary to use an external bypass for temperature control units without pressure limitation. In order to prevent the risk of overpressure in the system, the thermal fluid must always be adjusted to room temperature before switching off. This will prevent damage to the temperature control unit or the application. Any isolating valves must remain open (pressure equalization).

- The temperature and the dynamics of the process are determined by the flow temperature. A
 differential temperature (Delta T) forms between flow temperature and process temperature.
 This differential temperature may have to be limited, because Delta T might exceed permissible
 limits of the application (glass apparatus) and cause bursting. Adjust the Delta T value to your application.
- Select the thermal fluid to be used in such a way that it not only permits the minimum and maximum working temperature but is also suitable with regard to fire point, boiling point, and viscosity. In addition, the thermal fluid must be compatible with all the materials in your system.
- Avoid bending the temperature control and cooling water hoses (if required). Use suitable angle
 pieces and lay the hose connections with a large radius. Take the minimum bending radius from
 the data sheet of the temperature control hoses used.
- The selected hose connections must be able to withstand the thermal fluid, the working temperatures and the admissible maximum pressure.
- Check the hoses at regular intervals for any material fatigue (e.g. cracks, leaks).
- Keep the temperature control hoses as short as possible
 - The inside diameters of the temperature control hoses must correspond at least to the pump connections. Select bigger inside diameters for longer line lengths to compensate for pressure loss in the piping.
 - The viscosity of the thermal fluid determines the pressure drop and influences the temperature control result, especially at low working temperatures.
 - Too small connectors and couplers and valves can generate significant flow resistance.
 Your application will therefore be slower to reach its design temperature.
- Only valid for CC-520w, CC-525w, CC-906w and cooling baths used in continuous operation: Water as well as water and anti-freeze mixes must not be used as thermal fluids!
- Basically, you should only use the thermal fluid recommended by the manufacturer and only within the usable temperature and pressure range.



OPERATION MANUAL

- The application should be roughly at the same height of or below the temperature control unit if the temperature control is close to the boiling temperature of the thermal fluid.
- Fill the temperature control unit slowly, carefully and evenly. Wear the necessary personal protective equipment, such as goggles, heat-proof and chemically resistant gloves, etc.
- The temperature control circuit must be vented after filling and setting all required parameters. This is required to ensure trouble-free operation of the temperature control unit and hence your application.

INFORMATION

For water-cooled temperature control units, please take the cooling water temperature necessary for perfect operation and the required differential pressure from the data sheet. \rightarrow From page 94, section »Annex«.

"Pilot ONE®" controller 3.4

Please note the figure **"Pilot ONE"**. \rightarrow Page 6. The basic version of "Pilot ONE" (Basic) can be upgraded in three stages (from "Basic" to "Exclusive", from "Exclusive" to "Professional", and from "Professional" to "Explore").

Functional overview of "Pilot ONE®" 3.4.1

You can verify, and optionally upgrade, the delivery version of your temperature control unit using Pilot ONE, category "E-grade".

Overview of the E-grade variants	Temperature control units / E-grade	E-grade Basic	E-grade Exclusive	E-grade Profes- sional
	Brewing thermostats	-	х	0
	Unistat temperature control units	-	-	х
	UniCAL	-	-	х
	Other temperature control units	х	0	0
	 E-grade "Explore" (only for "Unistat"-series temperature control units) The E-grade includes E-grade "Professional" functionalities. In addition it includes: Sensor calibration for return sensor Return temperature display Display of performance (estimated or calculated from volume flow) Display mode "Explore" instead of "numerically large" PB commands correspond to DV-E-grade, except: Temperatures only in 0.01 °C and volume flows only in 0.1 l/min. 	_	_	0
	E-grade "DV-E-grade" – All interface commands are enabled. – Temperatures in 0.001 C, volume flows in 0.001 I/min possible (cf. E- grade Explore) Note: Only the interface commands are enabled, not the correspond- ing menu items in Pilot ONE!	0	0	0
	E-grade "OPC-UA" – OPC-UA interface via Ethernet. This E-grade additionally includes DV-E-grade functionalities.	0	0	0
	X = Standard equipment, O = optional, – = not possible			

Overview of E-grade functions

Function	E-grade Basic	E-grade Exclusive	E-grade Professi- onal
Temperature control			
Controller parameterization: predefined ¹ / TAC ²	X/-	-/X	-/X
Find parameters: Fast Identification / With Preliminary Test / Estimate Control Parameters	-/-/X	X/X/-	X/X/-
Sensor calibration for control sensor ³ : x -point	2	5	5
Monitoring: Low level and overtemperature ⁴	х	х	х
Adjustable alarm thresholds	х	х	х
VPC ⁵ (Variable Pressure Control)	х	х	х
Ventilation program	х	х	х
Automated compressor	х	х	х
Setpoint limitation	х	х	х
Program controller: x programs / max x steps	-/-	3/15	10/100
Ramp function: linear / non-linear	-/-	X/-	x/x
Temperature control mode: Internal and process	-	х	х
Adjustable max. heating and cooling capacity	-	х	х
Display & operation			
Temperature display: 5.7" Touchscreen	х	х	х
Display mode: graphically / numerically large / Explore	-/X/-	-/X/-	X/X/-
Display resolution: 0.1 °C / 0.01 °C	X/-	X/X	x/x
Graphic display for temperature curves: Window, full screen and scalable	х	х	х
Calendar, date and time	х	х	х
Language: CZ, DE, EN, ES, FR, IT, JP, KO, PL, PT, RU, TR, ZH	х	х	х
Temperature format switchable: °C, °F and K	х	х	х
Display mode (screen) can be switched by swiping	х	х	х
Favorites menu	х	х	х
User menu (Administrator Level)	-	-	х
2. Setpoint	-	-	х
Connections			
Digital interface RS232	х	х	х
USB interfaces: Host und Device	х	х	х
Ethernet RJ45 interface	х	х	х
Pt100 external sensor connection	_	х	х
External control signal (ECS STANDBY ⁶)	х	х	х

¹ TAC function available as a 30 day evaluation version.

² True Adaptive Control.

⁴ For temperature control units with integrated overtemperature protection.

- ⁵ For units with variable-speed pump or external bypass.
- ⁶ Standard on Unistats, otherwise via optional Com.G@te or POKO/ECS interface.

OPERATION MANUAL

³ Internal Pt100 and External Pt100.

OPERATION MANUAL

Chapter 3

Function	E-grade Basic	E-grade Exclusive	E-grade Professi- onal
Programmable potential-free contact (ALARM ¹)	х	х	х
AIF (Analog Interface) 0/4–20 mA or 0–10 V ²	X ³	х	х
Digital interface RS485 ⁴	х	х	х
Comfort & other			
Visual / audible alarm signal	х	х	х
Autostart (Power Failure Auto)	х	х	х
Plug & Play technology	х	х	х
Technology glossary	х	х	х
Remote control / data visualization via Spy software	х	х	х
E-grade evaluation versions available (valid for 30 days)	х	х	х
Storage/loading temperature control programs onto USB stick	-	х	х
Copying settings	_	-	х
Service data recorders (flight recorder)	х	х	х
PB commands ⁵	х	х	х
Communication Watchdog	_	-	х
Process data recording directly to a USB stick: Setpoint, Actual Value Internal and Actual Value Process / Heating Capacity %, Cooling Capaci- ty %, and Pump Pressure / Pump Speed, and VPC Pressure	-/-/-	X/X/-	x/x/x
Calendar start	_	-	х

3.5 Clock/event function

3.5.1 Rechargeable accumulator

"Pilot ONE" is fitted with a clock that continues to run even when the temperature control unit is switched off. The energy required for this purpose is provided by a rechargeable accumulator, which is automatically charged when the temperature control unit is switched on. The accumulator is dimensioned so that the clock can also continue to run for prolonged switch-off intervals (up to several months). If, after extremely prolonged switch-off time, time and date have been deleted, leaving the temperature control unit switched on for a few hours will usually suffice (no temperature control required). During this time, you can reset the time and date.

If after switching the unit off and back on again, the previously set time and date re-appear, it can be safely assumed that the rechargeable accumulator is defective. In this case, please contact Customer Support. \rightarrow Page 92, section **»Contact data«**.

3.5.2 Programmable event function

The "Calendar Start" offers a programmable event function. This enables you to enter a time at which the event is repeatedly triggered on a daily basis (until the activity in the menu is reset). 2 event types are currently selectable:

¹ Standard on Unistats, otherwise via optional Com.G@te or POKO/ECS interface.

² Via optional Com.G@te.

³ Restricted, see entries under "Temperature Control".

⁴ Via optional Com.G@te.

⁵ See Data Communication Manual. Anything controllable via the graphical user interface of the Pilot ONE can be controlled by PB commands.



3.5.2.1 Event function "Alarm clock event"

Several acoustic signals are used.

3.5.2.2 Event function "Program event"

After selecting **"Program event"** when configuring the event function, you will be prompted for the number of the program to be started. The program will be started automatically when the programmed event time is reached. If the temperature control unit is not active, this will also be started.

3.6 Operation via the touch screen


The entire operation is via the **>Touchscreen< [88]**. These functions can be activated by tapping the displayed text boxes/icons once. This also changes the display.

INFORMATION

You can cancel the current dialog or dialog sequence at any time by pressing the "ESC" touch button. When canceling a dialog or dialog sequence, it may be necessary to confirm the cancellation again. When canceling a dialog sequence, settings made earlier in the dialog sequence are discarded. Check your already carried out settings and re-enter as needed.

3.7 Display instruments

Display instruments

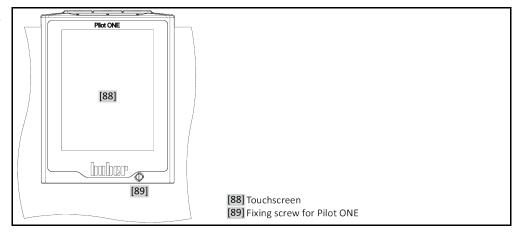

The following meters are available:


>Touchscreen< [88]

3.7.1 The touchscreen [88]

The most important display and operating instrument. Shows both standard variables (setpoint, actual value, setpoint thresholds...), and also menu guidance, error information output and operation.



3.8 Control instruments

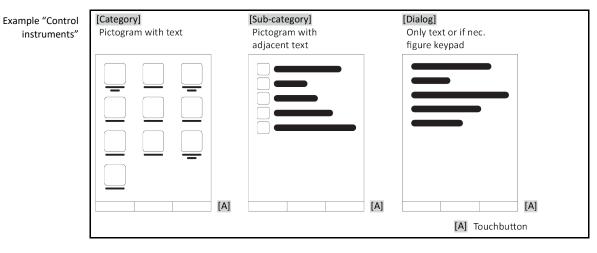

INFORMATION

To exit the "Categories Menu", sub-categories, menu items, press the "Home" touch button (house) or the arrow. After 2 minutes of inactivity, the category/ sub-category or the Favourites menu is automatically closed and you return to the "Home" screen. Dialogs are **not** canceled/closed after 2 minutes of inactivity.

3.8.1 The touchbuttons

Depending on the situation, the touch buttons can be assigned different functions. For example:

- Select the "Home" screen (house)
- Back (arrow to left)
- Favourites (star)
- Add to favorites (star with a plus sign)
- Select the "Categories menu" (menu)
- Confirm entry
- Start/stop
- etc.

3.8.2 The categories

For clarity we have grouped the Operation and Setting of Pilot ONE in various categories. A category is selected by tapping it.

3.8.3 The sub-categories

The sub-categories are parts of a category. This is where you will find the entries that we have grouped together for you in the selected category. Not all the categories also contain sub-categories. Tap on a sub-category to select it.

3.8.4 The dialogs

Tapping on a category or sub-category displays its dialogs. Dialogues may appear e.g. as text, a numeric or a alphanumeric keyboard. Dialogs allow you for example to enter settings or start created temperature control programs. Within a dialog, selection must always be confirmed with the "OK" touch button. When cancelling a dialog with the "ESC" touch button it may be necessary to confirm the cancellation again.



3.9 **Function examples**

3.9.1 Display of software version

PROCEDURE

- ➢ Go to the "Categories Menu".
- Tap on the category "System Settings". Tap on the category "System Info". \triangleright
- \geq
- > Tap on the sub-category "Software Version".

The software versions of the electronics will be displayed:

Display of software versions	Device name Operating temperature range	
	Serial number: xxxxx Load Code: xxxxxxxx Create Code: xxxxxxxxx	
	Pilot version: xxxxxxxxx.xx.xxxxxxxxx.x jjjj-mm-tt 00:00:00 Serial number: xxxxxxx	
	OS version: x.x	
	CAN controller: Vxx.xx.xxx	
	Bootloader version: xxxxxxxxx.xx.xx.xxxxxxx.x	
	Active Board: Vxx.xx.xx mmm tt jjjj 00:00:00 Serial number: xxxxxxx Active Board bootloader revision: x.x	

> Tap either the "ESC" or "OK" touchbutton. You will be returned to the previous screen.

> Tap on the "Home" touchbutton (house) to return to the "Home" screen.

3.9.2 Start & Stop

How to start and stop the temperature control process. Prerequisites: You have entered a setpoint.

PROCEDURE

➢ Go to the "Home" screen.

Start

- Tap on the "Start" touchbutton.
- Confirm the start of temperature control by tapping on "OK". The correct selection will be displayed graphically and temperature control will start immediately. If tapping on "OK" is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the "Home" screen again. Try to start the temperature control unit again.

Stop

- > Tap on the "Stop" touchbutton.
- Confirm the stop of temperature control by tapping on "OK". \triangleright The correct selection will be displayed graphically and temperature control will stop immediately and the pump with overrun for about 30 seconds. Wait until the pump stops. If tapping on "OK" is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the "Home" screen again. Try to stop the temperature control unit again.

INFORMATION

With the "Stop" touchbutton, you can also stop Degassing, Venting and Circulation. The requirement for this is that the corresponding task is active.



3.9.3 Copying the settings to a data carrier **Only valid for E-grade Professional**

The current settings of the temperature control unit are saved as a file to a connected USB flash drive. Use this file to restore the settings of a temperature control unit or to copy them to another temperature control unit. The models of the temperature control units used must be identical to copy the settings between temperature control units. You cannot transfer data between different models

You can also transfer the set value of the overtemperature protection to another temperature control unit. You must **always** verify this value on the respective temperature control unit and adjust it if necessary. Make sure that the value of the overtemperature protection has been adapted to the thermal fluid used.

3.9.3.1 Saving to a USB flash drive

PROCEDURE

- Insert a USB flash drive into the "USB 2.0 interface host". The USB flash drive must have at least 1 MB of free space.
- \triangleright Go to the "Categories Menu".
- Tap on the category "System Settings". Tap on the category "Copy Settings". \geq
- \triangleright
- Tap on the dialog entry "Save to USB flash drive".
- Confirm your choice by tapping on "OK". The content of the USB flash drive is displayed. If necessary, select the location (folder) on the USB flash drive.
- > Confirm the selection of the memory location by tapping on "OK".
- > Enter a name for the file to be saved. Confirm the entry by tapping on "OK". Alternatively, you can accept the proposed name by tapping on "OK".
- Read the message and confirm by tapping on "OK". The file with the settings was created on the USB flash drive.
- Remove the USB flash drive from the temperature control unit.

Loading from a USB flash drive 3.9.3.2

PROCEDURE

- Insert a USB flash drive with the saved file into the "USB 2.0 interface host".
- Go to the "Categories Menu".
- Tap on the category "System Settings".
 Tap on the category "Copy Settings".
- > Tap on the dialog entry "Load from USB flash drive".
- > Confirm your choice by tapping on "OK". The content of the USB flash drive is displayed.
- Select the file you want to load.
- > Confirm the selection of the file by tapping on "OK".
- \triangleright Select the setting group to be loaded from the list. A multiple selection is possible.
- \triangleright Confirm your choice by tapping on "OK".
- Read the message on the >Touchscreen< [88].
- Switch off the temperature control unit. The settings have been loaded to the temperature control unit.
- Remove the USB flash drive from the temperature control unit.

Restore factory settings 3.9.4

Use this function to reset the temperature control unit to various basic states. This is useful if you want to reverse various settings relatively quickly.

INFORMATION Restoring the factory settings is **only** possible if the temperature control unit is **not** carrying out a task. If a task is active, do not switch off the temperature control unit until the application allows this. Restoring the factory settings cannot be revoked. Depending on the type of resetting to factory settings performed you may have to enter parameters (process safety, thermal fluid used, cutout value etc.) again.

ററം



"X" = Value is reset, "-" = Value is not reset (A) = All together; (B) = Device parameters without OT; (C) = Device Parameter; (D) Program Control-ler; (E) = Menu; (F) = Com.G@te

Description	(A)	(B)	(C)	(D)	(E)	(F)
Acoustics		х	х	-	-	-
Display functions Warnings; Messages; Inactive Menu Items; Resolution Temperature	х	х	х	-	-	-
Display modes	х	х	х	-	-	-
AutoStart response	х	х	х	-	-	-
Limitations DeltaT limit; max. heating variable; max. cooling variable; max. permis- sible current consumption (for 230 V temperature control units)	x	x	x	_	_	_
Com.G@te (only if Com.G@te is connected) Analog Interface (configuration input)						
Meaning of AIF input; response on cable break; measurement range setting	х	х	х	-	-	х
Calibration	-	-	-	-	-	х
In case of analog error; current / voltage switching; filter constant	х	х	х	-	-	х
Configuration Output AIF output variable; Output range setting	х	х	х	-	-	х
Calibration	-	-	-	-	-	х
RS232/RS485 (Hardware RS, Baud rate, device address); ECS Standby; POKO Alarm	х	х	х	_	_	_
Properties (other) Change thermal fluid (thermal fluid, liter indicator, bypass usage); cooling bath / cooling (CC-E only)	x	х	x	_	-	-
Sensor calibration	-	-	-	-	-	-
Automated compressor	х	х	х	-	-	-
Enter program	х	-	-	х	-	-
Pump settings Pump speed setpoint; Pump pressure setpoint; Control Mode	х	х	х	-	-	-
Fluid properties Thermal fluid; liter indicator	х	х	х	-	-	-
Configuration manual parameters						
Protective function High Limit Alarm Internal; Low Limit Alarm Internal; High Limit Alarm Process; Low Limit Alarm Process; Hydrostatic Correction; Warning Time Limit (CC-E only); Min. Level (for temperature control units with analogue level sensor); Max. Level (for temperature controllers with analogue level sensor); Freeze Protection (optional)	x	x	x	_	_	_
Setpoint	х	х	х	-	-	-
Setpoint limitation Min. Setpoint; Max. Setpoint	х	х	х	-	_	-
Language	х	х	х	-	-	-
Temperature format	х	х	х	_	_	_
Temperature control mode	х	х	х	_	-	_
Favourites menu	х	-	-	-	х	-

Description	(A)	(B)	(C)	(D)	(E)	(F)
Pilot ONE Ethernet IP address; Subnet Mask; Remote Access	х	х	х	-	-	I
Time format	х	х	х	-	-	-
2. Setpoint	Х	Х	Х	-	-	-

3.9.4.1 Restore to factory settings without overtemperature protection

PROCEDURE

- ➢ Go to the "Categories Menu".
- Tap on the category "System Settings".
- > Tap on the category "Factory Settings".
- > Choose between the dialog entries "Unit Control Data without OT", "Menu", "Programmer" and "Com.G@te". These entries do **not** restore the overtemperature protection. Tap on the required dialog entry.
- Tap on the "OK" touchbutton to confirm your choice.
- Read the message displayed. Tapping on "Yes" restores the factory settings, tapping on "No" cancels the procedure. The message "Restart system!" appears on the >Touchscreen< [88].
- Switch the temperature control unit off. The selected control data have been reset.

3.9.4.2 Restore to factory settings including overtemperature protection

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "System Settings".
- Tap on the category "Factory Settings".
- Choose between the dialog entries "Unit Control Data" and "All together". These entries also reset the overtemperature protection. Tap on the required dialog entry.
- Tap on the "OK" touchbutton to confirm your choice.
- Read the message displayed. Tapping on "Yes" restores the factory settings, tapping on "No" cancels the procedure.

INFORMATION

In the following dialog, enter the overtemperature protection suitable to the thermal fluid used. To reset the cut-out values of the overtemperature protection to the factory setting from within the controller of Pilot ONE, enter for "Heater" 35 °C and for the expansion vessel 45 °C. "Process Safety" is factory set to "Stop" and is automatically reset to "Stop" when restoring the factory settings.

- Enter the temperature unit you want to use with Pilot ONE. The available choices are "Celsius" (°C)", "Kelvin (K)" and "Fahrenheit (°F)". Confirm your choice by tapping on "OK".
- \triangleright
- Read the message and confirm by tapping on "OK".
- Read the safety warning and confirm by tapping on "OK".
- Read the Note and confirm by tapping on "OK". \triangleright
- \triangleright Tap on the temperature format set by you in the controller (green text).
- \geq Enter the displayed red safety code using the numeric keypad that appears.
- Enter the value 35 °C for the "Heating OT value" using the numeric keypad that appears.
- \geq Confirm your entry by tapping on "OK".
- \triangleright Read the message and confirm by tapping on "OK".
- Read the safety warning and confirm by tapping on "OK". \triangleright
- Tap on the temperature format set by you in the controller (green text). \geq
- > Enter the displayed red safety code using the numeric keypad that appears.
- Enter the value 45 °C for the "Exp. vessel OT value" using the numeric keypad that appears.
- Confirm your entry by tapping on "OK". The message "Restart system!" appears on the \triangleright >Touchscreen< [88].
- Switch off the temperature control unit. The temperature control unit is reset.



4 Setup mode

4.1 Setup mode

CAUTION

Moving the temperature control unit during operation SERIOUS BURNS/FREEZING OF THE HOUSING PARTS/ESCAPING THERMAL FLUID Do not move temperature control units that are in operation.

4.1.1 Turning on the temperature control unit

PROCEDURE

➤ Turn on the temperature control unit using the >Mains switch< [37]. First, a system test is conducted to test the full functionality of the temperature control unit. All sensors, the all-pole disconnecting power disconnect relay for the power supply and the power electronics of the main heater and the main heater itself are checked. A message appears on the >Touchscreen< [88] in case of an error or if a warning is issued. In case of doubt, please contact Customer Support. → Page 92, section »Contact data«.</p>

INFORMATION The following data entry is required only for: a.) Initial start-up

- b.) resetting the temperature control unit to the factory settings. \rightarrow Page 48, section **»Restore** factory settings«.
- > After switching on the temperature control unit, tap on the required system language.
- ➤ Touch "OK" to confirm your selection.
- Tap on the thermal fluid used.
- Touch "OK" to confirm your selection.
- Read the message and confirm by tapping on "OK".
- > Enter the volume using the number keypad that appears.
- Touch "OK" to confirm your entry.
- Tap on the cooling bath used. This selection list is displayed only when the Pilot ONE was connected to a cooling bath.
- Touch "OK" to confirm your selection.
- Set the setpoint limits to the thermal fluid used. → Page 40, section »Information on the thermal fluids« and → Page 65, section »Setting the setpoint thresholds«.

INFORMATION If your temperature control unit is not integrated into a network, confirm the preset IP address (0.0.0.0) by tapping on "OK". This skips the network settings.


- > Enter the required IP address using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- > Enter the required Subnet Mask using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- > Tap on the required remote control mode.
- Confirm your choice by tapping on "OK".

4.1.2 Turning off the temperature control unit

PROCEDURE

- > Warm the thermofluid to room temperature.
- Stop the thermoregulation.
- Switch off the temperature control unit using the **>Mains switch<** [37].



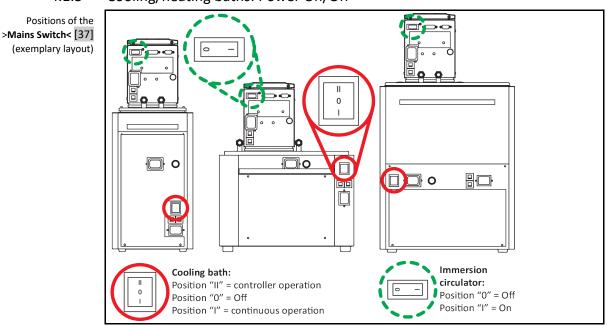

4.1.3 Cooling/heating baths: Power On/Off

4.1.3.1 Cooling/heating baths: Switching on (without an immersion circulator thermostat)

PROCEDURE

Switch on the cooling bath using the >mains switch< [37] (Position "I" continuous operation). The maximum cooling capacity of the cooling bath is permanently available in "continuous operation" (Position "I" of the >mains switch< [37]).</p>

4.1.3.2 Cooling/heating baths: Switching off (without an immersion circulator thermostat)

PROCEDURE

Switch off the cooling bath using the **>mains switch<** [37] (Position "0").

4.1.3.3 Cooling/heating baths: Switching on (with an immersion circulator thermostat)

INFORMATION The cold bath in combination with a hanger thermostat must only be used in switch position "II" - controller operation (on the cooling bath). The fuse may be tripped by excessive current consumption if ignored.

PROCEDURE

- Switch on the cooling bath using the >Mains switch< [37] (Position "II" controller operation).
- Switch on the Immersion circulator using the >mains switch< [37] (Position "I"). First, a system test is conducted to test the full functionality of the temperature control unit. All sensors, the all-pole disconnecting power disconnect relay for the power supply and the power electronics of the main heater and the main heater itself are checked. A message appears on the >Touchscreen< [88] in case of an error or if a warning is issued. In case of doubt, please contact Customer Support. -> Page 92, section »Contact data«.

INFORMATION The following data entry is required only for:

a.) Initial start-up
b.) resetting the temperature control unit to the factory settings. → Page 48, section »Restore factory settings«.

> After switching on the temperature control unit, tap on the required system language.

> Touch "OK" to confirm your selection.



- Tap on the thermal fluid used.
- Touch "OK" to confirm your selection.
- Read the message and confirm by tapping on "OK".
- > Enter the volume using the number keypad that appears.
- Touch "OK" to confirm your entry.
- Tap on the cooling bath used. This selection list is displayed only when the Pilot ONE was connected to a cooling bath.
- Touch "OK" to confirm your selection.
- Set the setpoint limits to the thermal fluid used. → Page 40, section »Information on the thermal fluids« and → Page 65, section »Setting the setpoint thresholds«.

INFORMATION If your temperature control unit is not integrated into a network, confirm the preset IP address (0.0.0.0) by tapping on "OK". This skips the network settings.

- > Enter the required IP address using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- > Enter the required Subnet Mask using the number keypad that appears.
- \succ Confirm your entry by tapping on "OK".
- > Tap on the required remote control mode.
- Confirm your choice by tapping on "OK".

4.1.3.4 Cooling/heating baths: Switching off (with an immersion circulator thermostat)

PROCEDURE

Switch off the hanger thermostat using the >mains switch< [37] (Position "0"). There is no need to switch off the cooling bath. Leave the >mains switch< [37] on Position "II" controller operation. If the temperature control unit is to be switched of for an extended period, set the >mains switch< [37] on the cooling bath to position "0" - Off.</p>

4.1.4 Setting the overtemperature protection

The overtemperature protection is not correctly set to the thermal fluid used MORTAL DANGER FROM FIRE

- > The overtemperature protection must be correctly set to the thermal fluid you are using.
- Always observe the safety data sheet of the thermal fluid.
- > Different working temperature ranges may apply for open and closed systems.
- > For Huber thermal fluids:
- For Huber thermal fluids, the maximum usable working temperature range is already indicated.
 Set the cut-out value of the overtemperature protection to the upper temperature limit of the
- working range of the thermal fluid.For thermal fluids of other manufacturers:
- Set the cut-out value of the overtemperature protection at least 25 K below the fire point of the thermal fluid.

INFORMATION

🚺 DANGER

The maximum settable cut-out value of the OT corresponds, for Huber thermal fluids, to the specified upper working temperature of the thermal fluid. The usable working temperature range may be less, if the overtemperature protection is correctly set. When controlling the temperature at the upper working temperature limit, the OT may be triggered due to tolerances.

4.1.4.1 General information on the overtemperature protection

The overtemperature protection is a device in the temperature control unit that works independently of the controller. The software and hardware is designed to test key functions and operating states during a self-test after the power supply has been switched on. If faults are detected, the enabling of the electrical assemblies in the temperature control unit is blocked. During operation, the sensors are tested for short-circuit and interruption.

The bath or flow temperature is monitored to ensure the safety of your system. It is set immediately after you have filled the system with thermal fluid.

Our temperature control units not only offer the possibility to specify the cut-out value of the overtemperature protection but also provide the possibility to specify the shutdown mode of the tem-

perature control unit. In a typical setting the temperature control unit switches off both the temperature control as well as the circulation (stop according to DIN EN 61010) after reaching the cut-out value. This monitors a possible defect in the control of the heater. A strong exothermic reaction near the cut-out value can also, under certain circumstances, result in a shutdown of the temperature control unit. In this case, however, a shutdown would be fatal.

Our temperature control units provide the possibility to work with the shutdown mode Process Safety. In this mode, temperature control (cooling) and circulation work as always. This is therefore a possibility to respond to exothermic reactions.

INFORMATION The default shutdown mode of the overtemperature protection is "Stop according to DIN EN 61010". When resetting to factory settings, the overtemperature protection is reset to the default shutdown mode "Stop according to DIN EN 61010"!

Upon delivery, the cut-out value of the overtemperature protection is set to 35 °C. An alarm is triggered by the temperature control unit shortly after turning on the power if the temperature of the thermal fluid just filled is higher than the cut-out value set for the overtemperature protection. Set the overtemperature protection to the thermal fluid you are using.

When setting a new cut-out value for the overtemperature protection, you will be prompted to enter a randomly generated and displayed code via the numeric keypad displayed. Only after successful entry will you be able to change the cut-out value.

4.1.4.2 Setting "OT limit: heating"

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Safety".
- Tap on the category "Overtemperature".
- Read the safety warning and confirm by tapping on "OK".
- Read the Note and confirm by tapping on "OK".
- \triangleright Tap on the dialog entry "OT limit: heating".
- \triangleright Confirm your choice by tapping on "OK".
- Tap on the temperature format set by you in the controller (green text).
- Enter the displayed red safety code using the number keypad that appears.
- Enter the value for the "Heating OT value" using the numeric keypad that appears. This value must be 25 K below the fire point of the thermal fluid you are using.
- Confirm your entry by tapping on "OK".

4.1.4.3 Setting "OT expansion vessel"

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Safety". Tap on the category "Overtemperature". \triangleright
- \triangleright Read the safety warning and confirm by tapping on "OK".
- ≻ Read the Note and confirm by tapping on "OK".
- > Tap on the dialog entry "OT limit: Expansion vessel".
- \triangleright Confirm your choice by tapping on "OK".
- \triangleright Tap on the temperature format set by you in the controller (green text).
- \triangleright Enter the displayed red safety code using the numeric keypad that appears.
- Enter the value for the "Exp. vessel OT value" using the number keypad that appears.
- Confirm your entry by tapping on "OK".

Setting "Process Safety" 4.1.4.4

You have two options.

"Stop" in accordance with DIN EN 61010

Once the overtemperature protection cut-out value is reached, the temperature control unit (heating, cooling circuit and circulation pump) switches off (default setting).



"Process Safety"

Once the overtemperature protection cut-out value is reached, the heater switches off and the cooling circuit and circulation pump continue to operate. In case of an emergency (possible exothermic reaction), the full cooling capacity is thus available. Please make sure that the automated compressor is set to Permanently ON ([System Settings] > [Power/ECO Settings] > [Compressor ON/OFF/AUTO] > [Permanently ON]).

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Safety". Tap on the category "Process Safety". \geq
- Read the safety warning and confirm by tapping on "OK".
- Choose between the modes "Stop" and "process safety".
- \triangleright Confirm your choice by tapping on "OK".

Monitoring via "Display OT values" 4.1.4.5

PROCEDURE

- Go to the "Categories Menu".
- \triangleright Tap on the category "Safety".
- Tap on the category "Display OT values". You will receive an overview of the current measured temperature value of the overtemperature protection sensor, the set cut-out values and the set cut-out mode (Process Safety). Some temperature control units have 2 overtemperature protection sensors, and so 2 values are shown for these temperature control units.
- Tap on the "OK" touchbutton after you have read/checked the information.

4.1.5 Testing overtemperature protection for functionality

DANGER Overtemperature protection (OT) does not trip MORTAL DANGER FROM FIRE Test the response of the device every month and after each change of the thermal fluid in \geq order to assure proper functioning. NOTE The steps below are carried out without permanent monitoring of the temperature control unit DAMAGE TO AND IN THE VICINITY OF THE TEMPERATURE CONTROL UNIT The following actions may only be carried out while constantly monitoring the temperature control unit and the application! INFORMATION Run the test only if the temperature of the thermal fluid used is about 20 °C. You must NOT leave the temperature control unit unattended during the entire test of the overtemperature protection. The Pilot ONE controller contains a description for carrying out the overtemperature protection test.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Safety".
- > Tap on the category "OT Test".
- The description for carrying out the overtemperature protection test is displayed.
- Tap on the "OK" touch button after you have read the information.

Adjusting the Delta T limiter 4.1.6





OPERATION MANUAL

INFORMATION The temperature dynamics inside the reactor/process temperature are dictated by the flow temperature. A differential temperature (Delta T) occurs between the flow temperature and the temperature inside the reactor. The greater the permitted Delta T, the better the energy transmission and hence the faster the speed reaches the setpoint. However, damage could result if the temperature difference limits are exceeded (bursting of the application e.g. glass apparatus). This difference in temperature may have to be restricted depending on the application (glass apparatus).

4.1.6.1 Changing the Delta T limiter

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Protection Options".
- Tap on the sub-category "Delta T limiter".
- Set the value of Delta T in line with the glass apparatus.
- Confirm your entry by tapping on the "OK" touchbutton.

4.2 The temperature control circuit


Every temperature control unit with a Pilot ONE control panel has its own PID controller for internal and process temperature control. For many temperature control tasks it is sufficient to use the factory-set control parameter. Our years of experience and current developments in control technology are used in these control parameters.

If a process control system is used, ideally send the specification of the temperature setpoint digitally to the temperature control unit. For this purpose, an Ethernet and USB port is available at the Pilot ONE and a RS232 interface at the temperature control unit. With the optional Com.G@te, an additional RS485 interface is added to your temperature control unit. Optionally you can integrate the temperature control unit in a Profibus environment. \rightarrow From page 73, section **»Interfaces and** software update«.

INFORMATION

The capacity adjustment of the temperature control unit is optimized so that the specified processes are run through in the fastest possible time. This increases the productivity of the system as a whole and thus saves energy in the long term.

Diagram of temperature control circuit

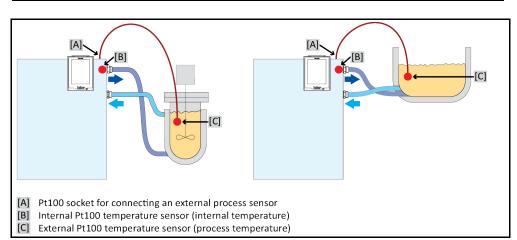

4.2.1 Select temperature control: Internal or process

PROCEDURE


- Go to the "Categories Menu".
- \geq
- Tap on the category "Temperature Control". Tap on the category "Process/Internal". Choose between the entries "Internal" and "Process (Cascade)". \geq
- Confirm your choice by tapping on "OK". \triangleright



4.2.2 Temperature control to internal temperature

With internal temperature control, a control circuit is used to control the temperature at the internal Pt100 temperature sensor. This Pt100 temperature sensor is built into the device and is located close to thermal fluid outlet (forward flow) or in the bath tank.

Representation of an optimum internal temperature control

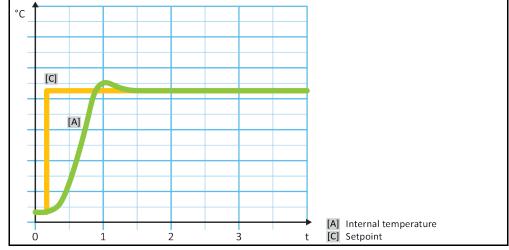

4.2.3 Temperature control to process temperature


Certain temperature control tasks require that the temperature is recorded elsewhere for best results than described. Setting the temperature to process temperature makes alternatives available. When setting the temperature to process temperature, an additionally connected external Pt100 temperature sensor is used in conjunction with a master controller (cascade controller). The internal sensor at the supply line is integrated with the slave controller. This temperature control method is used e. g. for the thermostatic control of jacket vessels. The setpoint setting is valid for the process controller. It in turn calculates a target value for the internal controller to optimally adjust the process setpoint.

NOTE

Incorrect installation of the process sensor (Pt100)

MATERIAL DAMAGE CAUSED BY FAULTY TEMPERATURE CONTROL

- > Disruption of measured value recording due to static build-up.
- The process sensor (Pt100) must have a screened supply cable.
- If the sensor tube is metallic, take care to avoid ground loops.
- The connection cable should not be unnecessarily long.
 - Make sure that the process sensor is properly attached at the measurement point and that there is good thermal coupling
- The sensor itself must have good insulation from the screen or the protective grounding (R > 20 MΩ).

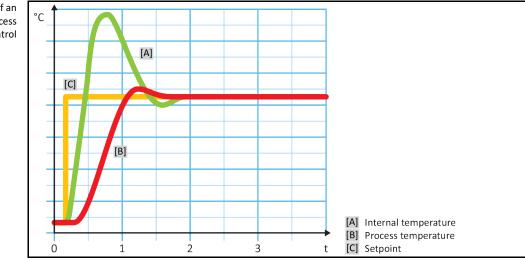

Representation of an optimum process temperature control \triangleright



OPERATION MANUAL

4.2.4 Delta T limiter

The Delta T limiter is a part of the temperature control which protects the system or the process. The Delta T limiter is given a limit value. The Delta T limiter then reacts when the limit value is reached during heating or cooling.

The temperature control mode "Process (Cascade)" evaluates the temperature difference between the flow temperature and the process temperature. The default setting of the limit value is 100 K. If the limit value and the temperature sensor are set properly, the load limits, e.g. from a glass apparatus, are not exceeded. As the limit value is approached, the cooling or heating capacity is adjusted. The DeltaT limiter is **not** a safety device.

4.2.5 Monitoring the Pt100 temperature sensors

The Pt100 temperature sensors are constantly monitored for their electrical status. If the status "Sensor faulty" occurs during temperature control, the temperature control process is stopped immediately and a device message is displayed. This applies for all temperature sensors connected within the temperature control unit.

4.2.6 Optimum control parameters for optimum temperature control

If the adjustment of the temperature does not correspond to the quality of the illustrations shown above, you can adjust the control parameters. With Huber temperature control units, there are various ways of finding the optimum control parameters. Depending on the facilities of the temperature control unit, you can choose the following processes:

- Use ex-factory parameters (standard) Estimate Control Parameters (only useful for bath thermostats with E-grade Basic and
 - internal control)
 - (from E-grade Exclusive)
- With Preliminary Test (from E-grade Exclusive)

4.2.7 Sub-category: "Select auto/expert mode".

Fast Identification

NOTE Use of the "Expert mode" without a thorough knowledge of I&C technology. MATERIAL DAMAGE TO THE APPLICATION > Only use this mode if you have a thorough knowledge of I&C technology.

> Here you can select whether the control parameters are set in the "Automatic mode" or in the "Expert mode". To enter settings in "Expert mode" requires in-depth knowledge of process measuring and control technology. Incorrect or insufficient settings can severely impair the function of the temperature control unit.

INFORMATION In "Expert mode", "Configuration auto" is deactivated and only "Configuration manual" is possible.

This is how to change mode:

PROCEDURE

- Go to the "Categories Menu".
- ≻
- Tap on the category "Temperature Control". Tap on the category "TAC/Manual". Tap on the sub-category "Select auto/expert". \geq
- Choose between the dialog entries "Automatic mode" and "Expert mode".
- Confirm your choice by tapping on "OK".

4.2.8 Sub-category: "Configuration auto"

INFORMATION You can only select this menu entry if "Automatic mode" is set.

4.2.8.1 Sub-category: "Find parameters"

4.2.8.1.1 Dialog entry: "Fast Identification"

Considering the small effort, the "Fast Identification" function of the controlled system provides you with very quick and reliable adapted control parameters. These control parameters achieve a fast and very accurate tuning performance. The more complex but also more precise identification "With Preliminary Test" is only very rarely required.

INFORMATION Do not make any changes to the temperature control unit and to the application (temperature control unit / external application) after your system has been started. Changes include e.g. filling / draining process chamber, changing the mixer speed, change of position of Pt100 process control sensor etc.

PROCEDURE

- Before setting the control parameters, make sure that the temperature control unit has reached the set setpoint and has been controlling the temperature at this setpoint for a few minutes. Do not stop the temperature control.
- Do not make any changes to the temperature control unit and the application during parameter determination.
- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- > Tap on the sub-category "Configuration auto".
- Tap on the sub-category "Find parameters".
- Tap on the dialog entry "Fast Identification".
- Confirm your choice by tapping on "OK".
- Read the message displayed and confirm it by tapping on "OK".
- Select the thermal fluid used from the list.
- Confirm your choice by tapping on "OK".
- Choose between the dialog entries "Internal" and "Process (Cascade)".
- Confirm your choice by tapping on "OK".
- Enter a new setpoint using the number keypad that appears. This should be at least 10 K away from the current setpoint.
- Confirm your entry by tapping on "OK". Setting of the control parameters using "Fast Identification" starts and after a while a message appears on the display.
- Read the message displayed and confirm it by tapping on "OK".

4.2.8.1.2 Dialog entry: "With Preliminary Test"

In some complex applications, the "Fast Identification" of the control circuit may not yet lead to optimum control. This can happen, in particular, if the hydraulic build-up cannot be adjusted to the necessary circulation volumes. \rightarrow Page 41, section **»To be noted when planning the test**«.

A further optimization of control behavior can be achieved if you select the control parameterization "With Preliminary Test". With this, the control parameters within the set limits of the minimum and maximum setpoint are determined. Temperature control is then also carried out to the setpoint limits in some circumstances.

INFORMATION

Before starting the automatic controller parameterization, verify the correct setting of the minimum and maximum setpoint. A limitation to the actually used operating temperaturerange is advantageous. Do not make any changes to the temperature control unit and to the application (temperature control unit / external application) after your system has been started. Changes include e.g. filling / draining process chamber, changing the mixerspeed, change of position of Pt100 process sensor etc.

Since the working temperature range can sometimes be very large, finding the parameters takes correspondingly longer in this mode. The control defines up to three temperature setpoints and automatically processes them one after another. One of these is below room temperature, one roughly at room temperature and one above room temperature, if the setpoint limits allow this.



PROCEDURE

- Before setting the control parameters, make sure that the temperature control unit has reached the set setpoint and has been controlling the temperature at this setpoint for a few minutes. Do not stop the temperature control.
- ⊳ Do not make any changes to the temperature control unit and the application during parameter determination.
- \geq Go to the "Categories Menu".
- \geq Tap on the category "Temperature Control".
- \geq Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration auto". Tap on the sub-category "Find parameters". \triangleright
- ⊳
- Tap on the dialog entry "With Preliminary Test".
- Confirm your choice by tapping on "OK". \triangleright
- Read the message displayed and confirm it by tapping on "OK".
- \triangleright Select the thermal fluid used from the list.
- \triangleright Confirm your choice by tapping on "OK".
- Choose between the dialog entries "Internal" and "Process (Cascade)". \geq
 - Confirm your choice by tapping on "OK". Setting the control parameters using "With Preliminary Test" starts and a message appears on the touchscreen after a while.
- Read the message displayed and confirm it by tapping on "OK".

4.2.8.1.3 Dialog entry: "Estimate Control Parameters"

Even with the simple temperature control units, we offer a further advantage in comparison with comparable bath thermostats available on the market. You can modify an existing control parameter by entering the thermal fluid used and the quantity of thermal fluid. This version is available with bath thermostats without connected external application.

INFORMATION

All the necessary technical data for the thermal fluids listed for Pilot ONE are stored in the controller. If the thermal fluid that you use is not shown in the list, choose the most similar thermal fluid with regard to temperature range and viscosity. \rightarrow Page 40, section »Information on the thermal fluids«.

PROCEDURE


- Before setting the control parameters, make sure that the temperature control unit has reached the set setpoint and has been controlling the temperature at this setpoint for a few minutes. Do not stop the temperature control.
- \triangleright Do not make any changes to the temperature control unit and the application during parameter determination.
- \geq Go to the "Categories Menu".
- Tap on the category "Temperature Control". Tap on the category "TAC/Manual".
- \triangleright
- Tap on the sub-category "Configuration auto". \geq
- Tap on the sub-category "Find parameters".
- Tap on the dialog entry "Estimate Control Parameters". \triangleright
- Confirm your choice by tapping on "OK". \triangleright
- \triangleright Read the message and confirm it by tapping on "Yes". The current control parameters will be overwritten.
- ≻ Select the thermal fluid used from the list.
- Confirm your choice by tapping on "OK". \triangleright
- \triangleright Read the message displayed and confirm it by tapping on "OK".
- Enter the filling volume using the number keypad that appears. \geq
- \geq Confirm your entry by tapping on "OK".
- Tap twice on the "Arrow" touchbutton to return to the category "Temperature Control". \triangleright
- Tap on the category "Process/Internal". \triangleright
- Choose between the dialog entries "Internal" and "Process (Cascade)". \triangleright
- \triangleright Tap on "OK" to confirm your choice.

4.2.8.2 Sub-category: "Control Dynamics"

You can choose between a faster tuning performance with a possible and accepted small temperature overshoot and a tuning performance without temperature overshoot. The default setting is "Fast, small overshoot".



The overshoot always refers to the leading temperature. If, for example, you have activated the process temperature control, this will be the leading temperature. Contrary to this, the bath or flow temperature must always be ahead of the process temperature. To achieve the best possible transmission of energy the largest possible temperature difference between bath and flow temperature and the process temperature is required. Note the illustration **"Representation of an optimum process temperature control"**. \rightarrow From page 57, section **"Temperature control to process temperature**". This can only ever be done with the largest possible volume flow of thermal fluid. In the setting **"Fast, small overshoot"**, the combination of high thermal fluid flow rate and the superbly designed control electronics hardly ever results in an overshoot, while reaching the setpoint as quickly as possible. Contrary to the mode **"Faster, small overshoot"** there is the setting **"Without overshoot"**. The target temperature is approached more cautiously and thus aperiodically. It takes longer to adjust to the selected setpoint. The statement "without overshoot" is only valid with little external disturbance influence. Please note the specifications. \rightarrow Page 41, section **"To be noted when planning the test"**.

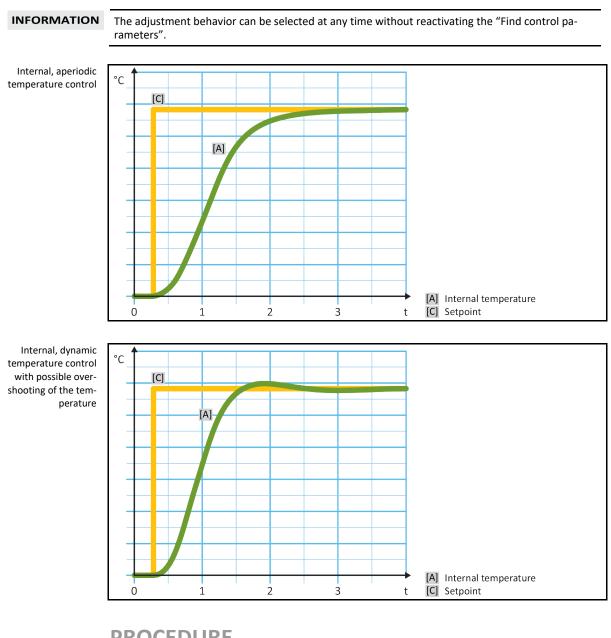

PROCEDURE

- ➢ Go to the "Categories Menu".
- > Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- > Tap on the sub-category "Configuration auto".

OPERATION MANUAL

Chapter 4

- Tap on the sub-category "Control Dynamics".
- ≻ Choose between the dialog entries "Fast, small overshoot" and "Without overshoot".
- \geq Confirm your choice by tapping on "OK".

4.2.8.3 Sub-category: "Fluid Properties"

4.2.8.3.1 Sub-category "Select Fluid"

Under this entry, you can select the thermal fluid used from a list.

PROCEDURE

- ➢ Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- > Tap on the category "TAC/Manual".
- Tap on the sub-category "Replandar".
 Tap on the sub-category "Configuration auto".
 Tap on the sub-category "Fluid Properties".
 Tap on the sub-category "Select Fluid".

- Select the thermal fluid used from the list.
- Confirm your choice by tapping on "OK".
- 4.2.8.3.2 Sub-category: "Bath/Circulation Volume"

Under this entry, you can enter the filling volume of the thermal fluid in your bath/circuit.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- > Tap on the category "TAC/Manual".

- Tap on the sub-category "Configuration auto".
 Tap on the sub-category "Fluid Properties".
 Tap on the sub-category "Bath/Circulation Volume".
- Read the message and confirm it by tapping on "OK".
- > Enter the filling volume using the number keypad that appears.
- > Confirm your entry by tapping on "OK".
- 4.2.8.3.3 Sub-category: "VPC/Bypass"

Under this entry, you can specify whether you use a bypass or not.

PROCEDURE

- ➢ Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration auto".
 Tap on the sub-category "Fluid Properties".
- Tap on the sub-category "VPC/Bypass".
- Choose between the dialog entries "Bypass used" and "Bypass not used".
- Confirm your choice by tapping on "OK".

4.2.8.3.4 Sub-category: "Show Fluid"

This entry provides you with an overview of the settings entered.

PROCEDURE

- ➢ Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- > Tap on the category "TAC/Manual".
- > Tap on the sub-category "Configuration auto".



- Tap on the sub-category "Fluid Properties". Tap on the sub-category "Show Fluid".
- \triangleright Tap on "OK" after you have read/checked the entries.

4.2.8.4 Sub-category: "Display parameters"

Here you can display the set parameters in "Automatic mode".

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
 Tap on the category "TAC/Manual".
- > Tap on the sub-category "Configuration auto".
- Tap on the sub-category "Display parameters".
- Tap on "OK" after you have read/checked the parameters.

Sub-category: "Configuration manual" 4.2.9

ΝΟΤΕ	Use of the "Expert mode" without a thorough knowledge of I&C technology.		
	MATERIAL DAMAGE TO THE APPLICATION		
	Only use this mode if you have a thorough knowledge of I&C technology.		
INFORMATION	In "Expert mode", "Configuration auto" is deactivated only "Configuration manual" is possible.		
INFORMATION	In "Expert mode", "Configuration auto" is deactivated only "Configuration manual" is possible.		

To enter settings in "Expert mode" requires in-depth knowledge of process measuring and control technology. Incorrect or insufficient settings can severely impair the function of the temperature control unit.

4.2.9.1 Sub-category: "Change parameters"

In this menu you manually configure the control parameters. If only the internal temperature is used as a target value, the control parameters are entered only under "Internal". If the process temperature is used as the target value, also the internal controller can be used, for example, when reaching the set point limit or in case of a Delta T limit. Consequently, parameter sets must be entered under all three items ("Internal", "Jacket" and "Process") when using the process temperature as a target value.

4.2.9.1.1 Sub-category: "Internal"

Enter the new values for "KP", "Tn" and "Tv" here one after the other.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration manual".
 Tap on the sub-category "Change parameters".
 Tap on the sub-category "Internal".

- > Enter the new "KP" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Enter the new "Tn" value using the number keypad that appears. \triangleright
- Confirm your entry by tapping on "OK". \geq
- Enter the new "Tv" value using the number keypad that appears.
- \triangleright Confirm your entry by tapping on "OK".

63

OPERATION MANUAL

Chapter 4

4.2.9.1.2 Sub-category: "Jacket"

Enter the new value for "KP" here.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control". Tap on the category "TAC/Manual".
- \geq
- > Tap on the sub-category "Configuration manual".
- Tap on the sub-category "Change parameters".
- Tap on the sub-category "Jacket".
- Enter the new "KP" value using the number keypad that appears. \geq
- Confirm your entry by tapping on "OK".
- Read the message and confirm it by tapping on "OK".

4.2.9.1.3 Sub-category: "Process"

Enter the new values for "KP", "Tn" and "Tv" here one after the other.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
 Tap on the category "TAC/Manual".
- > Tap on the sub-category "Configuration manual".
- Tap on the sub-category "Change parameters".
 Tap on the sub-category "Process".
- > Enter the new "KP" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- > Enter the new "Tn" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Enter the new "Tv" value using the number keypad that appears.
- \triangleright Confirm your entry by tapping on "OK".

4.2.9.2 Sub-category: "Display parameters"

In this function, the set manual parameters are displayed.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
 Tap on the category "TAC/Manual".
- > Tap on the sub-category "Configuration manual".
- Tap on the sub-category "Display parameters".
- \triangleright Tap on "OK" after you have read/checked the parameters.

4.2.9.3 Sub-category: "Control structure"

With this function, you have two different control structures available.

"Huber PID controller": Default setting

"Classic PID controller": This setting is exclusively used by Huber service engineers for service purposes.

PROCEDURE

- ➢ Go to the "Categories Menu".
- Tap on the category "Temperature Control".
 Tap on the category "TAC/Manual".
- > Tap on the sub-category "Configuration manual".



- Tap on the sub-category "control structure".
- Choose between the dialogue entries "Huber PID controller" and "Classic PID controller".
- Confirm your choice by tapping on "OK".

Sub-category: "Reset parameters" 4.2.10

With this function, you can reset the control parameters to the factory setting.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
 Tap on the category "TAC/Manual".
- Tap on the sub-category "Reset parameters".
- Read the message and confirm it by tapping on "Yes". The control parameters are reset/deleted. The temperature control unit can only be operated again after a complete restart.
- To do so, switch the temperature control unit off and back on again. The parameters have been reset.

Sub-category: "Display parameters" 4.2.11


In this function, the set parameters are displayed. Depending on the previous setting, these will be the "Automatic control parameters" or the "Manual control parameters".

PROCEDURE

- Go to the "Categories Menu".

- Tap on the category "Temperature Control". Tap on the category "TAC/Manual". Tap on the sub-category "Display parameters".
- Tap on "OK" after you have read/checked the parameters.

Setting the setpoint thresholds 4.2.12

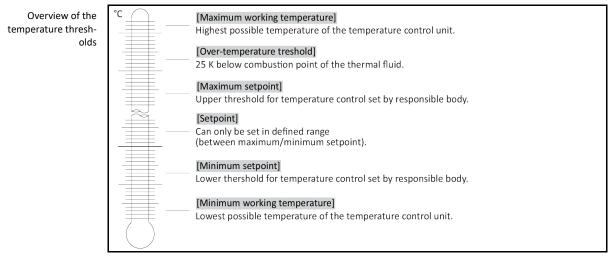

The limits for the minimum and maximum setpoint serve the safety of your system. They **must** be set for the application range of the thermal fluid before starting the first temperature control and when changing the thermal fluid. The maximum setpoint limit limits the setpoint setting for the bath or flow temperature. The minimum setpoint limit protects against high viscosity or freezing at low temperatures. The adjustable setpoint is then only available in the temperature range between the minimum and maximum setpoint limit.



PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Protection Options".
- Tap on the category "Setpoint Limits".
- Tap on the sub-category "Minimum Setpoint".
- > Enter the new value, using the numeric keypad that appears.
- Confirm your entry by tapping on "OK".
- In the display that follows, confirm your entry again by tapping on "OK". The correct selection will be displayed graphically and the "Min. setpoint" will be changed promptly. If tapping on "OK" is not correct, it will be displayed graphically for 2 seconds. The display will then return to the category "Setpoint Limits". Try changing the "Min. Setpoint" again.
- > Tap on the sub-category "Maximum Setpoint".
- > Enter the new value, using the numeric keypad that appears.
- Confirm your entry by tapping on "OK".
- In the display that follows, confirm your entry again by tapping on "OK". The correct selection will be displayed graphically and the "Max. Setpoint" will be changed promptly. If tapping on "OK" is not correct, it will be displayed graphically for 2 seconds. The display will then return to the category "Setpoint Limits". Try changing the "Max. Setpoint" again.

INFORMATION Check the set values for the minimum and maximum setpoint at any system change, especially when changing the thermal fluid.

4.2.13 Setting the setpoint

PROCEDURE

- Go to the "Home" screen.
- Tap on the keypad symbol next to "T_{Setpoint}".
- Enter a new setpoint using the number keypad that appears.
 The following must apply: [Minimum setpoint limit] ≤ [Setpoint] ≤ [Maximum setpoint limit].
 If these conditions are not met, a message will appear on the >Touchscreen< [88] and the entry will be ignored. In this case, delete the value just entered either with the "Arrow" key or with the "clear" key. Enter the setpoint again.
- Confirm your entry by tapping on "OK".
- In the display that follows, confirm your entry again by tapping on "OK". The correct selection will be displayed graphically and the setpoint will be changed immediately. If tapping on "OK" is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the "Home" screen again. Try changing the setpoint again.

4.3 Filling, venting, degassing and draining

Observe the wiring diagram. \rightarrow From page 94, section **»Annex«**.

 CAUTION
 Extremely hot / cold surfaces, connections and thermal fluids BURNS/FREEZING OF LIMBS
 Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
 Avoid direct contact with surfaces, connections and thermal fluids!
 Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).
 During an active circulation, the thermal fluid circuit is shut off by shut-off valves MATERIAL DAMAGE TO THE CIRCULATING PUMPS INSTALLED IN THE TEMPERATURE CONTROL UNIT
 Do not close the thermal fluid circuit during an active circulation by means of shut-off valves.
 Warm the thermal fluid to room temperature before stopping the circulation.



4.3.1 Filling, venting, degassing and draining the bath thermostat

CAUTION Non-compliance with the safety data sheet for the thermal fluid to be used

INJURIES

- Risk of injury to the eyes, skin, respiratory tract.
- The safety data sheet for the thermal fluid to be used must be read prior to using it and its content must be respected.
- > Observe the local regulations/work instructions.
- Wear your personal protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).
- Danger of slipping because floor and work area are contaminated. Clean the workplace; observe the proper disposal of thermal fluid and aids. → Page 17, section »Proper disposal of resources and consumables«.

4.3.1.1 Filling and venting the bath thermostat and the externally closed application

PROCEDURE

- Check that a hose was installed at the >Overflow< [12] (if present). The other end of the hose must be inserted in a suitable collecting container. Excess thermal fluid will leak at this point when the temperature control unit is overfilled. Hose and container must be compatible with the thermal fluid and its temperature.</p>
- ▶ Lift the **>Bath cover<** [93] from the temperature control unit.
- ➤ Using the filling accessories (funnel and/or beaker), carefully fill in suitable thermal fluid. While filling, pay attention to any necessary measures such as grounding the containers, funnels and other accessories. The thermal fluid flows from the temperature control unit via the hose connections to the external application. Follow the instructions for the proper disposal when cleaning filling accessories. → Page 17, section **»Proper disposal of resources and consumables**«.
- Switch on the temperature control unit.
- Go to the "Categories Menu".
- > Tap on the category "Temperature Control".
- Tap on the category "Start/Stop".
- Tap on the dialog entry "Start venting".
- Touch "OK" to confirm your selection.
- > Using the number keypad that appears, enter the time interval for venting. The default is 0.5 minutes.
- Confirm your entry by tapping on "OK" to start the filling process. The filling process is complete when the bath vessel is filled sufficiently. For temperature control units with heater, each heating coil of the heater must be below the thermal fluid level. This corresponds to the minimum level in the bath/minimal filling. There must be no air bubbles left in the external application (e.g. in the double jacket). The control electronics monitors the fill quantity and displays the level on the graphic display.
- Run the "venting program" for a few minutes. This ensures that trapped air bubbles can escape. These bubbles would result in a safety shutdown during an actual temperature control process.
- Stop venting. To do this, go to the category "Temperature Control".
- Tap on the category "Start/Stop".
- Tap on the dialog entry "Stop venting".
- > Touch "OK" to confirm your selection. Venting is stopped and the pump continues to run for approx. 30 seconds. Wait until the pump stops.
- Put the >Bath cover< [93] back onto the bath opening.

INFORMATION

The volume expansion of the thermal fluid depends on the working temperature range you wish to work in. Do not go below the minimum bath level/minimal level when working at the "lowest" working temperature and there should be no overflow from the expansion vessel/temperature control unit when working at the "highest" working temperature. In case of overfilling, drain the excess amount of thermal fluid. \rightarrow Page 69, section **»Draining the bath thermostat**«. Check if the thermal fluid can be reused. Follow the instructions for the proper disposal. \rightarrow Page 17, section **»Proper disposal of resources and consumables**«.

- Check the fill level in the collecting container. Empty the container when necessary and dispose of its contents properly.
- ➢ In case of overfilling, drain thermal fluid via the >Drain< 8 into a suitable container (see page → Page 69, section »Draining the bath thermostat«. Check if the thermal fluid can be reused. Follow the instructions for the proper disposal. → Page 17, section »Proper disposal of resources and consumables«.</p>

INFORMATION

ON Venting and degassing must be performed especially during commissioning and after changing the thermal fluid. This is the only way to ensure trouble-free operation. Follow venting with the degassing operation. → Page 68, section »Degassing of bath thermostat«.



4.3.1.2	Degassing of bath thermostat
CAUTION	 Hot or cold thermal fluid and surfaces BURNS TO LIMBS Avoid direct contact with the thermal fluids or the surfaces. Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).
FORMATION	When changing from low-boiling thermal fluid (low-boiling components) to higher boiling thermal fluids, remains of the low-boiling component may remain in the temperature control unit. Depending on the working temperature, the low-boiling component begins to boil and gas bubbles are formed that cause the pump pressure to momentarily collapse. This may cause a safety shutdown The gas bubbles reach the bath opening and can escape.
	If ice crystals form on the evaporator coil, water has accumulated in the thermal fluid. Degas if th is the case, to avoid damage to the temperature control unit.
	Thermal fluids are more or less hygroscopic (water-attracting). This effect increases, the lower the working temperature. The de-gassing mode below, which must be permanently monitored , also helps you remove any water residues from the temperature control circuit.
	PROCEDURE
	 Follow venting with the degassing operation. Prerequisite: You have filled and/or cleaned the temperature control unit in accordance with the instructions. → Page 67, section »Filling and venting the bath thermostat and the externally closed application« and/or → Page 85, section »Rinsing the thermal fluid circuit«. Go to the "Home" screen. Tap on the keypad symbol next to "Tsetpoint". Enter a setpoint using the number keypad that appears. This setpoint must be below the lower boiling thermal fluid. This setpoint will be increased in 10 K steps during the degassing process to the maximum working temperature. Touch "OK" to confirm your entry. In the display that follows, confirm your entry again by tapping on "OK". The correct selection will be displayed graphically and the "setpoint" will be changed immediate Iy. If tapping on "OK" is not correct, it will be displayed graphically for 2 seconds. The display we return to the "Home" screen. Try changing the setpoint again. Tap on the "Start" touch button. Confirm the start of temperature control by tapping on "OK". The correct selection will be displayed graphically and temperature control will start immediate Iy. If tapping on "OK" is not correct, it will be displayed graphically for 2 seconds. The display will return to the "Home" screen. Try starting the thermal control process again. Carry out temperature control to the entered setpoint until no more gas bubbles rise up. Increase the setpoint by 10 K and carry out temperature control until no more gas bubbles rise up. Repeat increasing the setpoint by 10 K until the maximum working temperature of the thermal fluid's maximum working temperature. Confirm the stop of the temperature control process by tapping on "OK". The correct selection will be displayed graphically. The temperature on the "Home" screen. Try stopping the temperature control process is stopped immediately and the pump continues to run fo

4.3.1.3 Draining the bath thermostat

SEVERE BURNS/FROSTBITE OF LIMBS

Hot or very cold thermal fluid

- > Before draining, ensure that the thermal fluid has room temperature (20 °C).
 - If, at this temperature, the thermal fluid is too viscous to be drained: Control the temperature of the thermal fluid for a few minutes until the viscosity will allow drainage. Never control the temperature of the thermal fluid when the drain is open.
 - Danger of burns when draining thermal fluids at temperatures above 20 °C.
- > Wear your Personal Protective Equipment for draining.
- Only drain using suitable drainage hose and collecting container. These must be compatible with the thermal fluid and its temperature.

PROCEDURE

Baths with >Drain valve< [4]

- Remove the knurled screw at the >Drain< [8].</p>
- Connect a suitable drain hose to the >Drain< [8].</p>
- Place the other end of the hose in a suitable container.
- > Open the >Drain valve< [4] by turning it counterclockwise (turn 90° left as far as it will go). The thermal fluid will flow from the external application via the bath vessel and the draining hose into the container. Check if the thermal fluid can be reused. Follow the instructions for the proper disposal. → Page 17, section »Proper disposal of resources and consumables«.</p>
- > Wait until the external application and the bath are empty.
- > Open the connection >Circulation flow< [1].
- > Open the connection >Circulation return< [2].
- Leave the temperature control unit open for a while to allow it to dry out and the residue to drain. Without screw caps and with the >Drain valve< [4] open.</p>
- Close the >Drain valve< [4] by turning it clockwise (turn 90° right as far as it will go).</p>
- Close the connection >Circulation flow< [1].</p>
- Close the connection >Circulation return< [2].</p>
- After drying out, remove the drain hose and re-fit the knurled screw to the >Drain< [8].</p>
- The bath is now drained.

Baths without >Drain valve< [4]

- Have a suitable container ready to catch the thermal fluid.
- ➢ Open the knurled screw at the >Drain< [8]. As soon as you have opened the knurled screw, the thermal fluid will flow from the external application over the bath and into the container. Wait until the external application and the bath are empty. Check if the thermal fluid can be reused. Follow the instructions for the proper disposal. → Page 17, section »Proper disposal of resources and consumables«.</p>
- > Open the connection >Circulation flow< [1].
- > Open the connection >Circulation return< [2].
- Leave the temperature control unit open for a while for the residue to fully drain and to allow it to dry out (without screw caps).
- Close the connection >Circulation flow< [1].</p>
- Close the connection >Circulation return< [2].</p>
- Re-fit the knurled screw to the >Drain< [8].</p>
- The bath is now drained.

5 Normal operation

5.1 Automatic operation

Extremely hot / cold surfaces, connections and thermal fluids BURNS/FREEZING OF LIMBS

- Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
- Avoid direct contact with surfaces, connections and thermal fluids!
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).

NOTE

NOTE

During an active circulation, the thermal fluid circuit is shut off by shut-off valves MATERIAL DAMAGE TO THE CIRCULATING PUMPS INSTALLED IN THE TEMPERATURE CONTROL UNIT

- > Do not close the thermal fluid circuit during an active circulation by means of shut-off valves.
- Warm the thermal fluid to room temperature before stopping the circulation.

5.1.1 Temperature control

5.1.1.1 Starting the temperature control process

The temperature control process can be started after filling and complete venting.

PROCEDURE

- Go to the "Home" screen.
- Tap on the "Start" touchbutton.
- Confirm the start of temperature control by tapping on "OK". The correct selection will be displayed graphically and temperature control will start immediately. If tapping on "OK" is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the "Home" screen again. Try to start the temperature control unit again.

5.1.1.2 Ending the temperature control process

When the temperature control unit is switched off, the thermal fluid temperature is higher/lower than room temperature

DAMAGE TO THE TEMPERATURE CONTROL UNIT AND THE GLASS APPARATUS/APPLICATION

- > Bring the thermal fluid up to room temperature using the temperature control unit.
- > Do not close the shut-off valves in the thermal fluid circuit.

The temperature control can be ended at any time, the pump continues to run for approx. 30 seconds. The compressor is switched off after the stepper motor valve has reached a defined position to control the cold output.

PROCEDURE

- Go to the "Home" screen.
- > Tap on the "Stop" touchbutton.
- Confirm the stop of temperature control by tapping on "OK".
 - The correct selection will be displayed graphically and temperature control will stop immediately and the pump with overrun for about 30 seconds. Wait until the pump stops. If tapping on "OK" is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the "Home" screen again. Try to stop the temperature control unit again.

INFORMATION

The compressor is not switched off until the stepper motor valve has reached a defined position. The status line [Field 10] displays the relevant information.



5.1.2 Temperature control via a created temperature control program

5.1.2.1 Starting the temperature control program

A temperature control program can be started after filling and complete venting.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Programmer/Ramp".
- Tap on the category "Start/stop Program".
- > Tap on the dialog entry of the temperature control program to be started.
- Confirm your choice by tapping on "OK".
- Read the message and confirm it. You temperature control unit will start the temperature control program and the temperature control programmed in it will start.
- Read the Note and confirm by tapping on "OK".

5.1.2.2 Ending/cancelling the temperature control program

NOTE

When the temperature control unit is switched off, the thermal fluid temperature is higher/lower than room temperature

DAMAGE TO THE TEMPERATURE CONTROL UNIT AND THE GLASS APPARATUS/APPLICATION

- Bring the thermal fluid up to room temperature using the temperature control unit.
- > Do not close the shut-off valves in the thermal fluid circuit.

You can either end the temperature control automatically using the parameters predefined in the temperature control program or end/interrupt the temperature control at any time manually. The temperature control is switched off immediately after, the pump continues to run for approx. 30 seconds. The compressor is switched off after the stepper motor valve has reached a defined position to control the cold output.

Manual ending/cancelling

PROCEDURE

- Go to the "Home" screen.
- Tap on the "Stop" touchbutton.
- Confirm the stop of temperature control by tapping on "OK".

The correct selection will be displayed graphically and temperature control will stop immediately and the pump with overrun for about 30 seconds. Wait until the pump stops. If tapping on "OK" is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the "Home" screen again. Try to stop the temperature control unit again.

INFORMATION

The compressor is not switched off until the stepper motor valve has reached a defined position. The status line [Field 10] displays the relevant information.

5.2 Cooling/heating baths: Automatic operation (without an immersion circulator thermostat)

Extremely hot / cold surfaces, connections and thermal fluids

BURNS/FREEZING OF LIMBS

- Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
- Avoid direct contact with surfaces, connections and thermal fluids!
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).


NOTE

During an active circulation, the thermal fluid circuit is shut off by shut-off valves MATERIAL DAMAGE TO THE CIRCULATING PUMPS INSTALLED IN THE TEMPERATURE CONTROL UNIT

Do not close the thermal fluid circuit during an active circulation by means of shut-off valves.
 Warm the thermal fluid to room temperature before stopping the circulation.

5.2.1 Cooling/heating baths: Temperature control

Positions of the >Mains Switch< [37] (exemplary layout)

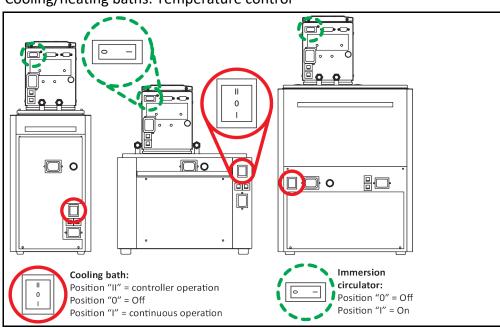

5.2.1.1 Cooling/heating baths: Starting the temperature control process

The temperature control process can be started after filling and complete venting.

PROCEDURE

Turn on the cooling bath (Position "I" - continuous operation). The maximum cooling capacity of the cooling bath is permanently available in "continuous operation" (Position "I").

5.2.1.2 Cooling/heating baths: Ending the temperature control process

NOTE

When the temperature control unit is switched off, the thermal fluid temperature is higher/lower than room temperature

DAMAGE TO THE TEMPERATURE CONTROL UNIT AND THE GLASS APPARATUS/APPLICATION

Bring the thermal fluid up to room temperature using the temperature control unit.
 Do not close the shut-off valves in the thermal fluid circuit.

The temperature control can be ended at any time, the pump continues to run for approx. 30 seconds.

PROCEDURE

■ ➤ Turn off the cooling bath (Position "0"- continuous operation).



6	Interfaces and software update
ΝΟΤΕ	The specifications of the interface used are not being met. PROPERTY DAMAGE
	 Only connect components that meet the specifications of the interface used.
INFORMATION	The use of PB commands is described in our "Data communications PB" manual. This manual can be downloaded from www.huber-online.com.

6.1 Interfaces at the "Pilot ONE®" controller

NOTE

The Pilot ONE controller is not operated behind a firewall

PROPERTY DAMAGE

- Operate the controller Pilot ONE exclusively behind a firewall, if the local subnet is connected to the Internet or another risk-prone network.
- The best available technology is to be applied in order to provide sufficient security for the LAN!

Standard interfaces at the "Pilot ONE" top side

6.1.1

10/100 Mbps Ethernet for RJ45 network sockets

This is a fast and flexible interface. Standard 10/100 Mbps interface (Fast Ethernet), can be connected to any existing Ethernet network. Because this interface can also be connected to very large networks, the IT "Best Practices" (firewall) must be observed.

Usage:

Also - to be able to communicate with the "Pilot ONE" controller - the communication enable must be issued. This is an additional safety feature that prevents persons - possibly unintentionally - connecting to the wrong machine and implementing incorrect temperature control specifications. The following restrictions are possible:

- Deactivated
- Always on (PLC)
- 12h Inactivity Timer
- 10min Inactivity Timer

If, for example, "10min Inactivity Timer" is selected, the connection must be made within 10 minutes after confirmation at the control. If this does not happen, the connection is refused.

INFORMATION Communication with the Pilot ONE is via TCP (Transmission Control Protocol), Port 8101. The interfaces used must meet the specifications of the generally accepted standards.

6.1.2 USB-2.0 interface

INFORMATION The interfaces used must meet the specifications of the generally accepted standards.



OPERATION MANUAL

Chapter 6

6.1.2.1 USB-2.0 interface, host


USB-2.0 connection (for connector A), e.g. for data memories.

6.1.2.2 USB-2.0 interface, device

USB-2.0 connection (for Mini-B connector) for communicating with a computer.

6.2 Cooling/heating baths: Interfaces at the back

Interfaces at the back of the cooling bath (exemplary layout)

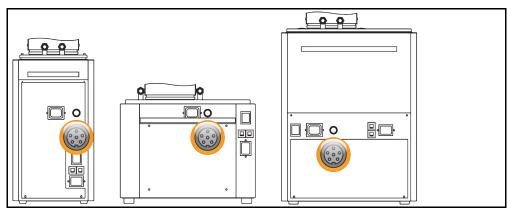


6.2.1 Activation connector



This connector is used to connect the cooling bath with the hanger thermostat. This enables control of the cooling bath via the hanger thermostat.

6.3 Cooling/heating baths: Interfaces on the hanger thermostat

Interfaces at the back of the immersion circulator (exemplary layout)

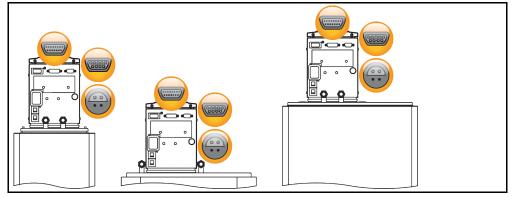

6.3.1 Service interface

This interface is exclusively used by Huber service engineers for service purposes. An adapter cable makes this interface a RS232 serial port.

6.3.2 RS232 (option) / activation of cooling bath

This connector is used to connect the hanger thermostat with the cooling bath. This enables control of the cooling bath via the hanger thermostat.

INFORMATION

The interfaces used must meet the specifications of the generally accepted standards.

6.3.3 Connection jack for Pt100 process control sensor

A temperature sensor located in the connected application (Pt100, 4-wire technology, Lemosa connector) is connected to the Pt100 connection jack. The external actual temperature is then recorded and the operating temperature of the temperature control unit is permanently calculated and adjusted.


OPERATION MANUAL

INFORMATION

Depending on the operating temperature, isolation losses and exothermic heat, the operating temperature (flow temperature) of the application can be significantly less than the setpoint of the application. In this context, safety-critical thresholds for the temperature control fluid must be strictly observed.

The control results contained in the data sheet can only be achieved with **shielded** sensor leads. We recommend the external Pt100 process control sensor from the Huber accessories program.

Pin assignment (front view)

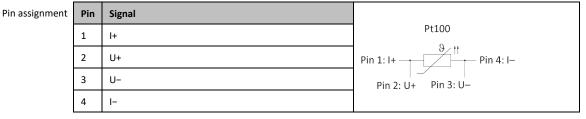

6.4 Interfaces at the Com.G@te[®] (optional)

NOTE

Connecting to the interfaces at the temperature control unit during operation DAMAGE TO THE INTERFACES

- When devices in operation are connected with interfaces of the temperature control unit, interfaces may get damaged.
- Before connecting, ensure the temperature control unit and the device to be connected are turned off.

Com.G@te external interfaces: LEVEL, POKO, AIF, ECS,RS232/RS485 (left to right).









If the **>Com.G@te**< [46] is not connected, proceed as follows.

PROCEDURE

- Switch off the temperature control unit.
- Plug the >Com.G@te< [46] into the slot provided.</p>
- Switch on the temperature control unit. "Unistat Control ONE" detects the new component. The >Com.G@te< [46] is operable.</p>

INFORMATION Note that you may only change >**Com.G@te**< [46] when the temperature control unit is turned off. In the category "Interfaces" you can change the settings of each function such as PoKo, analog interface and RS232/RS485.

You may order **>Com.G@te<** [46] with Order No. 6915. You can extend an existing **>Com.G@te<** [46] via a connecting cable (Order No. 16160).

6.4.1 Jack LEVEL (Com.G@te[®] external only)

For level monitoring in the **>Sight glass<** [23].




This connection enables you to connect an external float switch (Order No. 6152), which is positioned in the **>Sight glass<** [23], for monitoring the level of your externally closed application. Activation via a potential-free contact.

INFORMATION

The interface is specified as a digital input. Do not apply voltage or current.

Pin assignment (front view)

Pi

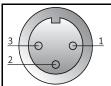

Pin assignment	Pin Description	
	1	Level test (bypass via contact $2 \rightarrow$ "absence")
2 Level – (GND)		Level – (GND)
	3	Level + (normally open contact)

6.4.2 Connector POKO (floating contact) alarm

Signal contact for external monitoring.



Notice the functions provided by PoKo in the category "Interfaces". The potential-free contact (PoKo) signals the state of the temperature control unit via the contact position. A closed operating contact means readiness for operation. If a fault or an error occurs, the operating contact is opened (this applies to the make contact between pin 1 and pin 2).

The following settings are offered:

- **"Off":** The POKO displays the OK state whenever the temperature control unit is ready for operation. The unit is ready for operation about 30 seconds after switching on, once the internal controller check has been completed. The OK status is terminated by switching off the unit or in case a fault occurs.
- "Internal temperature relative": Set the upper and lower limit of a temperature range around the setpoint by using the POKO "min. value" and the POKO "max. value".
 The potential-free contact indicates the condition, that the setpoint differs from the actual value by a smaller difference than was defined by the temperature range. Exceeding the set range only results in the potentialfree contact to switch (away from the OK state), but does not result in further reactions of the temperature control unit. The contact returns to the OK status if the actual value is within that range.

- "External alarm": The POKO relay becomes active (OK state) in this function only, if the temperature control unit is switched on and is in "Fault" mode. The advantage is that the alarm will not be triggered if the temperature control unit is switched off. Please use the POKO function "OFF" if the alarm function is wanted together with the operating current principle.
- "Unipump/PCS": This POKO function and the enable contact of the pump guarantee that the external pump runs synchronously with the pump in the temperature control unit if an external pressure booster pump is used in your temperature control circuit / cooling water circuit, i.e. the POKO assumes the OK state as soon as the internal pump is started up.

PCS: The POKO is used to communicate the state of the temperature control unit to the process control system.

POKO condition **ON** means the pump is running.


POKO condition **OFF** means the pump is not running, the temperature control unit is in stand-by mode.

- "External control": The PB commands "vPoKoExtMode" and "vPoKoState" use an external interface (Ethernet, RS232, TS485, USB device) to switch the POKO on and off. Please also refer to the software offered by us and the Data Communication Manual.
- "Process temperature relative": Set the upper and lower limit of a temperature range around the setpoint by using the POKO "min. value" and the POKO "max. value".
 The potential-free contact indicates the condition, that the setpoint differs from the actual value by a smaller difference than was defined by the temperature range. Exceeding the set range only results in the potential-free contact to switch (away from the OK state), but does not result in further reactions of the temperature control unit. The contact returns to the OK status if the actual value is within that range.
- "Unipump with Echo": This function is used to check whether the Unipump controlled by the POKO runs synchronous with the Unistat pump. For this purpose, the operating mode of the Unipump is applied via a make contact to the "LEVEL" socket. A fault is generated if asynchronous. This mode is useful if a Unipump is to be monitored either to guarantee the desired circulation or to prevent an unintended heating of the thermal fluid.
- **"Programmer":** This activates the individual segments associated during program creation with the POKO switching states.
- "Internal temperature absolute": This enables you to set a temperature range relative to the internal sensor (absolute temperature). The POCO is active outside of this range; the POCO is in-active inside of this range.
- "Process temperature absolute": This enables you to set a temperature range relative to the process temperature (absolute temperature). The POCO is active outside of this range; the POCO is inactive inside of this range.
- "Solenoid valve flow / return flow": This function is used to control a connected solenoid valve. It takes 60 seconds before the POKO switches on after you have started the pump in the temperature control unit. For example, the POKO is switched off before the pump has been completely stopped when the temperature control process / circulation is stopped. Thus, the POKO is switched on only as long as the full pump pressure exists.
- "Cooling": This function is used to open the cooling water supply with a solenoid valve only when the temperature control unit requires cooling water. The POKO is switched on when the temperature control unit is cooling.
- "Message available": The POKO switches on when a message is available at the Pilot ONE. This message can be a fault, warning or general message.
- **"Automatic filling":** The POKO switches on when the level has dropped below the minimum value. The POKO is switched off after 20 seconds when the level has risen above the minimum value. The POKO switches off immediately when the maximum value has been reached.

INFORMATION With th

With the floating contact, use only sheathed lines! The interface is specified as a digital output.

Pin assignment (front view)

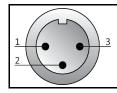

The connection is designed as a potential-free changeover contact. Closing contact between pin 1 and pin 2. Opening contact between pin 2 and pin 3. Contact load: 1 A at 24 V DC.



OPERATION MANUAL

Chapter 6


6.4.3 Jack AIF Reg-E-Prog



The analog interface has a programmable input channel and 3 output channels.

The analog interface of Com.G@te® is programmed in the category "Interfaces".

Pin assignment (front

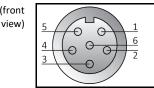

Pin assignment

t	Pin	Description	Signal	
	1	Current output, T external	0 ^{a)} /4 - 20 mA or 0 - 10 V	If using 0 - 10 V, integrate a
	2	Current output, setpoint	0 ^{a)} /4 - 20 mA or 0 - 10 V	500 Ω resistor.
	3 GND for analog outputs		GND	
	4	Analog input (programmable)	0 ^{a)} /4 - 20 mA or 0 - 10 V	Power input: 200 Ω working resistance Voltage input: 100 kΩ input resistance
	5	Current output, freely programmable	0 ^{a)} /4 - 20 mA or 0 - 10 V	If using 0 - 10 V, integrate a 500 Ω resistor.
	6 GND for analog input		GND	
^{a)} Consult with our Customer Support.				

6.4.4 Jack ECS (External Control Signal) standby

Enable signal ECS (external control signal) for starting/stopping the temperature control process.



Activation via a potential-free contact. Contacts 1 and 3 are internally bypassed. **ECS** is energized when E1 and E2 are connected by an external floating contact. Contact specification: min. 0.1 A / 24 V DC.

The functionality of the **ECS** is determined via the "Interfaces" category.

The following variants are offered:

- "No Action": Switching the contacts open/closed or closed/open has no effect.
- "Switching to second setpoint": A change from a closed to an open contact replaces the set setpoint with the value of the second setpoint. This altered setpoint is not fixed to the second setpoint but can be changed by the operator at the temperature control unit at any time. A change from an open to a closed contact does not cause any change and the temperature control process is not reset to the original setpoint.
- **"2. setpoint selective":** An open contact causes a thermoregulation to its original setpoint. A closed contact causes a thermoregulation to the second setpoint.
- "Internal / Process": If the contact is open, the device-internal temperature sensor is used for control. If the contact is closed, the additionally connected external temperature sensor is used for control. Switching the contacts open/closed or closed/open switches between these two temperature sensors.
- **"Temperature control on/off"**: The temperature control process starts when switching from an open to a closed contact. The temperature control process switches off when switching from a closed to an open contact.
- "Release": The temperature control process switches off, if the temperature control process is
 active and the contact is switched from closed to open. The temperature control process is not
 switched on, if contact is further switched from open to closed!
- "Reset Messages": All messages in the Pilot ONE are reset (if possible) if the enable signal changes from "open" to "closed". Resetting is delayed in the event of a fault until the temperature con-

trol unit is in stand-by mode. Messages can be reset as often as required, but faults can be reset only 3 times.

- "Program 1 Start/Stop": The temperature control program 1 is started if the enable signal changes from "open" to "closed". The temperature control program is stopped when changing from "closed" to "open".
- "Fill level alarm": When switching from closed to open contact, an alarm is triggered if the fill level is too low. For this purpose, a level sensor must be connected and placed in the > sight glass< [23] or in the external application.
- "Fill level warning": When switching from closed to open contact, a warning is generated if the fill level is too low. For this purpose, a level sensor must be connected and placed in the > sight glass< [23] or in the external application.

INFORMATION	The interface is specified as a digital input. Do not apply voltage or current.
-------------	---

Pin assignment (front view)

Pin assignment	Pin	Signal
	1.3	E2
	2	E1

6.4.5 Jack RS232/RS485 serial



A PC, a SPS or a Process Control System (PCS) can be connected to this jack for remote control of the controller electronics. Alternatively, a connection to a RS485 bus is possible. Before plugging in the cable, check the settings in the "Interfaces" category and adjust if necessary.

INFORMATION

The interfaces used must meet the specifications of the generally accepted standards.

Pin assignment (front view)



1

6

Pin assignment Pin

t	Pin	Signal	Description			
	Wirir	Wiring RS232				
	2	RxD	Receive Data			
	3	TxD	Transmit Data			
	5	GND	Signal GND			
	Wiring RS485					
	6	A with 120- Ω terminating resistor	-			
	7	A	-			
	8	В	-			

6.5 Firmware update

An instruction for running a firmware update can be found at www.huber-online.com.

Q

7 Service/maintenance

7.1 Messages from the temperature control unit


Messages output by the temperature control unit can be divided into various classes.

Follow the instructions displayed on the **>Touchscreen**< [88]. Once a message has been acknowledged, a symbol is output on the **>Touchscreen**< [88]. Tapping the symbol takes you to an overview of all messages in chronological order.

Displayed symbol:

7.2 Electrical fuse

OPERATION MANUAL

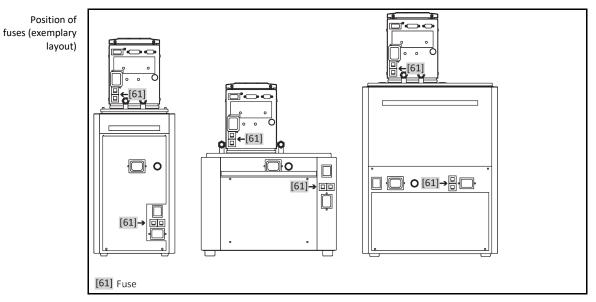

The thermal overcurrent circuit breakers for all-pole disconnection (L and N) are located at the back. In case of a fault (no function and/or no display) please first check if the overcurrent circuit breakers have tripped. If the overcurrent circuit breakers trigger again immediately after reversing, please unplug the power cord and contact Customer Support immediately. \rightarrow Page 92, section **»Contact data«**.


7.3 Replacing the "Pilot ONE®"

Replacing the electronics while the temperature control unit is operating MORTAL DANGER FROM FIRE

- Stop an ongoing temperature control process.
- Turn off the temperature control unit.
- > Also disconnect the temperature control unit from the power supply.

In case of fault, you can replace the "Pilot ONE" yourself. If you have questions or difficulties, please contact your dealer, your local office or our Customer Support.



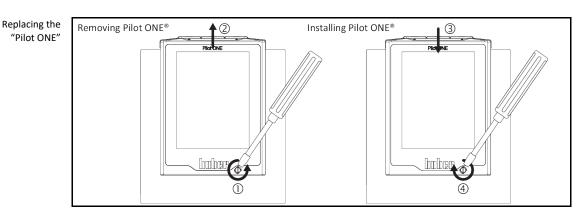

PROCEDURE

- > Turn off the temperature control unit.
- > Disconnect the temperature control unit from the power supply.
- > Release the >Fixing screw for Pilot ONE< [89] at the front of the housing.
- Carefully pull the "Pilot ONE" away upwards.
- ➤ Carefully insert the replacement "Pilot ONE".
- > Close the **>Fixing screw for Pilot ONE<** [89] at the front of the housing.
- > Connect the temperature control unit to the power supply.
- > Turn on the temperature control unit.

7.4 Maintenance

Cleaning/maintenance while the temperature control unit is operating MORTAL DANGER FROM ELECTRIC SHOCK

- > Stop an ongoing temperature control process.
- Turn off the temperature control unit.
- Also disconnect the temperature control unit from the power supply.

NOTE

N

\Lambda DANGER

Performing maintenance work not described in these operation manual

- MATERIAL DAMAGE ON THE TEMPERATURE CONTROL UNIT
- Please contact Huber for maintenance work that is not described in these operation manual.
- Maintenance work not described in these operation manual is reserved for qualified specialists trained by Huber.
- Safety-relevant components may only be replaced by equivalent ones. The specified safety values for the respective component must be observed.

7.4.1 Function check and visual inspection

Aonitoring intervals Cooling* Description		Maintenance interval	Comment	Person re- sponsible	
	A/W	Visually inspect hoses and hose connections	Prior to switching on the tempera- ture control unit	Exchange leaking hoses and hose connections prior to switching on the temperature control unit. → Page 82, section »Replacing temperature control or coolant hoses«.	Operating company and / or operators
	A/W	Check the fill level in the collecting container at the > Overflow < [12] (if present)	Prior to switching on the tempera- ture control unit	Check the fill level in the collecting container. Empty as required. Follow the instructions for the proper disposal of thermal fluid. → Page 17, section »Proper dis- posal of resources and consuma- bles «.	Operating company and / or operators
	A/W	Inspection in ac-	In accordance	→ Page 20, section »Temperature	Operating



Cooling*	Description	Maintenance interval	Comment	Person re- sponsible	
	cordance with the F-Gas Directive	with the F-Gas Directive	control units with fluorinated greenhouse gases/refrigerants«.	company	
A/W	Check the power supply cable	Prior to switching on the tempera- ture control unit or on relocation	Do not start the temperature control unit if the power supply cable is damaged.	Qualified electrician (BGV A3)	
A	Clean the perforat- ed sheet	As required	Clean the perforated sheet of the temperature control unit with a damp cloth	Operating company	
A/W	Thermal fluid inspection	As required	-	Operating company and / or operators	
A/W	Check the mechan- ical seals	Monthly	→ Page 87, section »Inspect the mechanical seal«.	Operating company and / or operators	
A	Check the liquefier fins	As required, after 3 months at the latest	→ Page 83, section »Clean liquefi- er fins (air-cooled temperature control unit)«.	Operating company and / or operators	
W	Check the hat-type strainer (dirt trap)	As required, after 3 months at the latest	→ Page 84, section »Clean hat- type strainer (dirt trap) (water- cooled temperature control unit)«.	Operating company and / or operators	
A/W	Overtemperature protection (OT) - functional check	Every month or after changing the thermal fluid.	→ Page 55, section »Testing overtemperature protection for functionality«.	Operating company and / or operators	
A/W	Check the tempera- ture control unit for damage and stabil- ity	Every 12 months or after a change of location	-	Operating company and / or operators	
W	Check the cooling water quality	Every 12 months	Descale the cooling water circuit as required. Documentation on water quality is available at: www.huber-online.com	Operating company and / or operators	
A/W	Exchange safety- relevant electric and electrome- chanical compo- nents	20 years	Have the exchange only carried out by certified personnel (such as Huber service engineers). Please contact Customer Support. → Page 92, section »Contact data« .	Operating company	
*A = Air cooling; W = Water cooling; U = Applicable only for Unistats					

7.4.2 Replacing temperature control or coolant hoses

Replace defective temperature control and/or coolant hoses **before** turning on the temperature control unit.

7.4.2.1 Replacing temperature control hoses

PROCEDURE

- ➢ Drain the temperature control unit and accessories. → Page 69, section »Draining the bath thermostat«.
- Replace defective temperature control hoses. Follow the instructions for the proper disposal.
 Page 17, section »Proper disposal of resources and consumables«.
- > Reconnect your external application. → Page 35, section »Connecting externally closed application«.
 > Fill the temperature control unit with thermal fluid. → Page 67, section »Filling and venting the bath thermostat and the externally closed application«.



- ➤ Vent the temperature control unit. → Page 67, section »Filling and venting the bath thermostat and the externally closed application«.
- Restart the temperature control unit in normal mode.

7.4.2.2 Replacing coolant hoses

PROCEDURE

- > Drain the cooling water. \rightarrow Page 90, section **»Draining the cooling water**«.
- ➢ Replace the defective coolant hoses. Follow the instructions for the proper disposal. → Page 17, section »Proper disposal of resources and consumables«.
- Reconnect the temperature control unit to the building's cooling water supply. → Page 31, section »Temperature control units with water cooling«.
- Restart the temperature control unit in normal mode.

7.4.3 Clean liquefier fins (air-cooled temperature control unit)

CAUTION Manual cleaning

RISK OF BEING CUT ON THE LIQUEFIER FINS

- → Wear suitable cut-resistant gloves for cleaning work.
- Depending on the ambient conditions, use cleaning equipment such as vacuum cleaners and/or a hand brush/brush. Follow the local regulations when cleaning. Do not clean the liquefier fins in a clean room with items like a brush and do not use a vacuum cleaner without an extra-fine particle filter.

Cleaning using pointed or sharp-edged tools
DAMAGE TO THE LIQUEFIER FINS

Clean the liquefier fins using suitable cleaning appliances.

INFORMATION

NOTE

Make sure there is adequate ventilation (removal of waste heat, fresh air supply) for the temperature control unit, in case of **air cooling, maintain wall clearance**. \rightarrow Page 23, section **»Exemplary illustrations of the cooling variants**« and \rightarrow Page 27, section **»Ambient conditions**«. The liquefier fins must be cleaned (dust) from time to time as only then will the temperature control unit perform at its maximum cooling capacity.

Identify the position of the ventilation grille, usually located on the front. With some temperature control units, the ventilation grilles on the side wall, rear or on the underside (table units) of the temperature control unit.

PROCEDURE

Ventilation grille on the front/rear or on a side wall

- Turn off the temperature control unit.
- Disconnect the temperature control unit from the power supply.
- > Remove the ventilation grille to create unhindered access to the liquefier fins.
- Clean the liquefier fins using suitable cleaning appliances. Observe the local regulations and ambient conditions when selecting cleaning appliances.
- Make sure the liquefier fins are not damaged or deformed as this will impair the air flow.
- Re-mount the ventilation grille after cleaning work.
- Connect the temperature control unit to the power supply.
- Turn on the temperature control unit.

PROCEDURE

Ventilation grille on the underside (table-top units)

NOTE

Cleaning of liquefier fins at the underside when the temperature control unit is filled DAMAGE CAUSED BY THERMAL FLUID PENETRATING THE TEMPERATURE CONTROL UNIT

Empty the temperature control unit before cleaning the liquefier fins at the underside of the temperature control unit.

OPERATION MANUAL

NOTE

- Switch off the temperature control unit.
- Disconnect the temperature control unit from the power supply.
- ➢ Drain the thermal fluid from the temperature control unit. → Page 69, section »Draining the bath thermostat«.
- ▶ Tilt the temperature control unit to remove the grille (if available) in front of the liquefier fins.
- Clean the liquefier fins using suitable cleaning appliances. Observe the local regulations and ambient conditions when selecting cleaning appliances.
- Make sure the liquefier fins are not damaged or deformed as this will impair the air flow.
- Re-mount the ventilation grille after cleaning work.
- Connect the temperature control unit to the power supply.
- ▶ Refill the temperature control unit with thermal fluid. \rightarrow Page 67, section »Filling and venting the bath thermostat and the externally closed application«.

7.4.4 Clean hat-type strainer (dirt trap) (water-cooled temperature control unit)

Building's isolating valves are not closed

PROPERTY DAMAGE CAUSED BY FLOODING OF ROOMS

- > Close the building's isolating valves in the cooling water supply and return lines.
- For table-top models, place a collection container below the >Cooling water drain< [15]. Observe the wiring diagram: → From page 94, section »Annex«.</p>

INFORMATION The strainer at the cooling water inlet must be inspected and cleaned on a regular basis, depending on water quality.

PROCEDURE

Table-top models:

- > Turn off the temperature control unit.
- Disconnect the temperature control unit from the power supply.
- Close the customer's shut-off valves in the cooling water supply and return lines.
- Place a collecting container below the >Cooling water inlet< [13].</p>
- Remove the cooling water supply line and take out the hat-type strainer for inspection and cleaning.
- Clean the suction strainer under running water.
- Following inspection/cleaning, reinsert the hat-type strainer and fasten the cooling water supply line.
- > Remove the collecting container from below the **>Cooling water inlet<** [13].
- Open the customer's shut-off valves in the cooling water supply and return lines.
- > Connect the temperature control unit to the power supply.
- Turn on the temperature control unit.

PROCEDURE

Free-standing models:

- Switch off the temperature control unit.
- Disconnect the temperature control unit from the power supply.
- > Close the customer's isolating valves in the cooling water supply and return lines.
- Remove the paneling around the cooling water supply [13], [14] and [15] (if present).
- Place a collecting container below the >Cooling water inlet< [13] and another collecting container below the >Cooling water drain< [15] (if present).</p>
- Open the ball valve of the >Cooling water drain< [15] (if present). If the temperature control unit is not equipped with a >Cooling water drain< [15]: Open the >Cooling water outlet< [13]. The cooling water will begin to drain out. Allow the cooling water to fully drain.
- Disconnect the >Cooling water inlet< [13] from the building's cooling water supply. Located immediately behind the >Cooling water inlet< [13] is the dirt trap.</p>
- Carefully detach the cover (hexagonal).
- Remove the metal strainer located below.
- Clean the metal strainer under running water.
- Re-insert the metal strainer after cleaning work.
- Carefully fasten the cover (hexagonal).
- Connect the >Cooling water inlet< [13] to the building's cooling water supply.</p>
- Close the ball valve of the >Cooling water drain< [15] (if present).</p>

Service/mair	itenance
Chapter 7	OPERATION MANUA
	 Remove the collecting containers below the >Cooling water inlet< [13] and below the >Cooling water drain< [15] (if present). Dispose of the contents of the collecting containers. Follow the instructions for the proper disposal. → Page 17, section »Proper disposal of resources and consumables«. Mount the paneling around the cooling water supply [13], [14] and [15] (if present). Open the customer's isolating valves in the cooling water supply and return lines. Connect the temperature control unit to the power supply. Switch on the temperature control unit.
INFORMATION	We are also happy to offer you service training. Please contact Customer Support. \rightarrow Page 92, section »Contact data« .
7.5	Thermal fluid inspection, replacement and circuit cleaning
	Observe the wiring diagram. \rightarrow From page 94, section »Annex« .
	Extremely hot / cold surfaces, connections and thermal fluids
	 BURNS/FREEZING OF LIMBS Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
	 Avoid direct contact with surfaces, connections and thermal fluids! Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).
ΝΟΤΕ	During an active circulation, the thermal fluid circuit is shut off by shut-off valves
	 MATERIAL DAMAGE TO THE CIRCULATING PUMPS INSTALLED IN THE TEMPERATURE CONTROL UNIT Do not close the thermal fluid circuit during an active circulation by means of shut-off valves. Warm the thermal fluid to room temperature before stopping the circulation.
7.5.1	Thermal fluid inspection
	Thermal fluid is not inspected on a regular basis
	 BURNS DUE TO REDUCED BOILING POINT Regularly check your thermal fluid whether it meets the specifications in the safety data sheet.
ΝΟΤΕ	Thermal fluid is not inspected on a regular basis
	DAMAGE TO THE HEAT EXCHANGER AND/OR ELECTROMECHANICAL PARTS.
	Regularly check your thermal fluid whether it meets the specifications in the safety data sheet.
INFORMATION	Oxidation Oxidation ages the thermal fluid and changes its characteristics (e.g. reduced boiling point). The reduced boiling point may cause overflow of very hot thermal fluids at the >Expansion vessel< [18] during temperature control of high temperatures. It may cause serious burns of the limbs. Hygroscopy During continuous temperature control below room temperature the thermal fluid accumulates water by hygroscopicity over time. Such a liquid mixture causes the evaporator to burst during temperature control in the minus range. This is caused by the water in the liquid mixture, which forms ice crystals on the evaporator. A temperature control of such a liquid mixture at high tem- peratures reduces the boiling point. A temperature control at high temperatures may cause an overflow of very hot thermal fluid at the >Expansion vessel< [18] as a result of the reduced boiling point. It may cause serious burns of the limbs.
7 5 0	Dincing the thermal fluid circuit

7.5.2 Rinsing the thermal fluid circuit

DANGER

Setpoint and overtemperature protection are not adjusted to the thermofluid MORTAL DANGER FROM FIRE

- > The cut-out value of the overtemperature protection **must** be adapted to the thermofluid. Set the
- cut-out value of the overtemperature protection 25 K below the fire point of the thermofluid.
 The setpoint set during rinsing **must** be adjusted to the thermofluid used.



UPERATION MANUA	L Chapter 7
	 Non-compliance with the safety data sheet for the thermal fluid to be used INJURIES Risk of injury to the eyes, skin, respiratory tract. The safety data sheet for the thermal fluid to be used must be read prior to using it and its content must be respected. Observe the local regulations/work instructions. Wear your personal protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear). Danger of slipping because floor and work area are contaminated. Clean the workplace; observe the proper disposal of thermal fluid and aids. → Page 17, section »Proper disposal of resources and consumables«.
ΝΟΤΕ	 Mixing different thermofluids in a thermal fluid circuit PROPERTY DAMAGE ➢ Do not mix different types of thermofluid (such as mineral oil, silicone oil, synthetic oil, water, etc.) in a thermofluid circuit. ➢ The thermal fluid circuit must be rinsed when changing from one type of thermal fluid to another. No residues of the previous type of thermal fluid may remain in the thermal fluid circuit.
Example: Connecting a short circuit hose	The inner components of the temperature control unit must be dried out Need to avoid boiling retardation during future uses (e.g. use of a silicone oil at temperatures above about 100 °C).
INFORMATION	 PROCEDURE > Drain the temperature control unit. → Page 69, section »Draining the bath thermostat«. Residual thermal fluid can remain in the pump chamber and the internal lines after draining. Leave
	 the temperature control unit with open valves for a while. Leave the drain hose mounted to the >Drain< [8]. Check the fill level in the collecting container at the other end of the drain hose. Follow the instructions for the proper disposal of thermal fluid. → Page 17, section »Proper disposal of resources and consumables«. Close the drain valves on the temperature control unit by turning them clockwise (turn 90° right as far as they will go). Connect the >Circulation flow< [1] with the >Circulation return< [2] at the temperature control using a bypass hose.
INFORMATION	 Perform the following steps without attaching a short circuit hose, if the application used by you (externally closed) is also dirty. In this case, leave your externally closed application connected to the temperature control unit. This rinses the temperature control unit and your application at the same time. Fill the system (minimum fill level) with the thermal fluid you wish to use. → Page 67, section »Filling and venting the bath thermostat and the externally closed application«. Vent the system. → Page 67, section »Filling and venting the bath thermostat and the externally closed application«.
	 Adjust the setpoint and the cut-out value of the overtemperature protection to the thermal fluid used. → Page 66, section »Setting the setpoint« and → Page 53, section »Setting the overtemperature protection«. Go to the "Categories Menu". Tap on the category "Temperature Control". Tap on the category "Start/Stop". Tap on the dialog entry "Start temperature control process". Touch "OK" to confirm your selection. The length of rinsing depends on the level of soiling. Tap on the category "Start/Stop".

Service/mail	ntenance	huber
Chapter 7		OPERATION MANUAL
	 ➤ Touch "OK" to confirm your selection. The temperature control proc > Drain the temperature control unit. → Page 69, section »Draining th > Repeat the steps "Filling", "Venting", "Start/Stop temperature control until the drained thermal fluid remains clear. > Remove the bypass hose after completely draining the temperature 	e bath thermostat«. ol process" and "Draining"
INFORMATION	Leave an application connected, if you have simultaneously rinsed a use closed).	ed application (externally
	 > Leave the >Drain< [8] open for a while to allow the thermal fluid rem control unit to evaporate. > Close the >Drain< [8] once the residual thermal fluid has evaporated > Dismount the drain hose. > Remove the collecting container. > Dispose of the collecting container, including its contents. Follow the disposal. → Page 17, section »Proper disposal of resources and conse > Re-connect your application. (Only if you have rinsed the thermal fluid hose.) 	instructions for the proper sumables«.

- \succ Fill the temperature control unit with thermal fluid. \rightarrow Page 67, section **»Filling and venting the** bath thermostat and the externally closed application«.
- ▶ Vent the temperature control unit. → Page 67, section **»Filling and venting the bath thermostat** and the externally closed application«. An externally open application does not need to be vented.
- Start the function "Degassing". → Page 68, section **»Degassing of bath thermostat**«. An externally open application does not need to be de-gassed.
- Restart the temperature control unit in normal mode.

7.6 **Cleaning the surfaces**

Extremely hot / cold surfaces, connections and thermal fluids

- **BURNS/FREEZING OF LIMBS**
- > Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
- Avoid direct contact with surfaces, connections and thermal fluids!
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).

NOTE

Exposed plug contacts

DAMAGE CAUSED BY FLUID INGRESS

- Protect unused plug contacts with the protective caps supplied.
- Clean surfaces only with a damp cloth.

A standard stainless steel cleaning agent is suitable for cleaning the stainless steel surfaces. Carefully clean painted surfaces (damp only) using a solution of sensitive-fabrics detergent. Observe the proper disposal of thermal fluid and aids. \rightarrow Page 17, section »Proper disposal of resources and consumables«.

7.7 Inspect the mechanical seal

NOTE

No visual inspection of the mechanical seal

MATERIAL DAMAGE IN THE TEMPERATURE CONTROL UNIT CAUSED BY LEAKING MECHANICAL SEAL

- Check the mechanical seal once a month.
- If case of leakage, stop the temperature control unit and contact Customer Support. \rightarrow Page 92, section »Contact data«.

Expect the formation of drops at the mechanical seal when operating with thermal fluids that evaporate only very slowly, as mechanical seals are never absolutely tight. These drops must be removed if necessary. \rightarrow Page 81, section **»Function check and visual inspection«**. The tightness of the mechanical seal must be visually checked. In case of a leakage, more thermal fluid escapes at the bottom of the temperature control unit. Observe the proper disposal of thermal fluid. \rightarrow Page 17, section »Proper disposal of resources and consumables«.



OPERATION MANUAL

Chapter 7

7.8 Plug contacts

NOTE

Exposed plug contacts DAMAGE CAUSED BY FLUID INGRESS

- > Protect unused plug contacts with the protective caps supplied.
- Clean surfaces only with a damp cloth.

Protective caps are supplied for all plug contacts. Make sure that any plug contacts not required are protective with the caps.

7.9 Decontamination/repairs

CAUTION

Returning a not decontaminated temperature control unit for repair

PHYSICAL INJURY AND PROPERTY DAMAGE CAUSED BY HAZARDOUS MATERIALS IN OR ON THE TEMPERATURE CONTROL UNIT

- > Carry out appropriate decontamination.
- > The decontamination process depends on the type and quantity of the materials used.
- Consult the relevant safety data sheet.
- > You will find a prepared return receipt at www.huber-online.com.

As the responsible body you are responsible for carrying out decontamination **before** third-party personnel come into contact with the temperature control unit / accessory. Decontamination must be carried out **before** the temperature control unit / accessory is returned for repair or inspection. Attach a clearly visible written notice stating that the temperature control unit / accessory has been decontaminated.

To simplify the process, we have prepared a form for you. This is available for download at www.huber-online.com.



8 Shutting down

8.1	Safety instructions and basic principles
A DANGER	 Connection/adjustment to the power supply not carried out by an electrician and/or connection to a power socket without protective earth (PE) MORTAL DANGER FROM ELECTRIC SHOCK Have the connection/adjustment to the power supply carried out by an electrician. Always connect the temperature control unit to safety sockets (PE).
DANGER	 Damaged power cable/power cable connection MORTAL DANGER FROM ELECTRIC SHOCK > Do not start up the temperature control unit. > Isolate the temperature control unit from the power supply. > Have the power supply cable/power supply connection replaced and inspected by an electrician. > Do not use a power cable that is longer than 3 m.
	Risk of tipping due to unstable temperature control unitSERIOUS INJURY AND PROPERTY DAMAGE➤ Avoid risk of tipping due to unstable temperature control unit.
CAUTION	 Non-compliance with the safety data sheet for the thermal fluid to be used INJURIES Risk of injury to the eyes, skin, respiratory tract. The safety data sheet for the thermal fluid to be used must be read prior to using it and its content must be respected. Observe the local regulations/work instructions. Wear your personal protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear). Danger of slipping because floor and work area are contaminated. Clean the workplace; observe the proper disposal of thermal fluid and aids. → Page 17, section »Proper disposal of resources and consumables«.
CAUTION	 Hot or very cold thermal fluid SEVERE BURNS/FROSTBITE OF LIMBS Before draining, ensure that the thermal fluid has room temperature (20 °C). If, at this temperature, the thermal fluid is too viscous to be drained: Control the temperature of the thermal fluid for a few minutes until the viscosity will allow drainage. Never control the temperature of the thermal fluid when the drain is open. Danger of burns when draining thermal fluids at temperatures above 20 °C. Wear your Personal Protective Equipment for draining. Only drain using suitable drainage hose and collecting container. These must be compatible with the thermal fluid and its temperature.
INFORMATION	All safety instructions are important and must be followed accordingly during working operations!

8.2 Switch-off

PROCEDURE

- > Turn off the temperature control unit.
- Disconnect the temperature control unit from the power supply connection.

8.3 Draining the temperature control unit

PROCEDURE

➢ Drain the temperature control unit. → From page 66, section »Filling, venting, degassing and draining«.

8.4 Draining the cooling water

INFORMATION	This section must be observed when using water-cooled temperature control units.
-------------	--

8.4.1 Draining process

Pressurized cooling water connections

RISK OF INJURY

- > Wear your personnel protective equipment (e.g. safety goggles).
- Carefully open the cooling water connection. Open slowly (1-2 signal edges) and drain the cooling water slowly.

NOTE

The building's isolating valves are not closed DAMAGE BY ROOM FLOODING

- Close the building's isolating valves in the cooling water supply and return lines.
 - For table-top models, place a collection container below the **>Cooling water outlet**< [14] and/or **>Cooling water drain**< [15] (if any).

PROCEDURE

Temperature control units with >Cooling water drain< [15]

- Close the isolating valves in the water supply on the temperature control unit (if present) and on the building side.
- Place a collecting container below the cooling water outlet [13], [14] and [15].
- Open the >Cooling water drain< [15] and disconnect the >Cooling water outlet< [14] from the water return. The cooling water will begin to drain out. Allow all the cooling water to drain out to prevent the risk of freezing during transport and storage!</p>
- Disconnect the >Cooling water inlet< [13] from the water supply.</p>
- Close the >Cooling water drain< [15].</p>

PROCEDURE

Temperature control units without >Cooling water drain< [15]

- Close the isolating valves in the water supply on the temperature control unit (if present) and on the building side.
- Place a collecting container below the cooling water outlet [13] and [14].
- Disconnect the >Cooling water outlet< [14] from the water return flow. The cooling water will begin to drain out. Allow all the cooling water to drain out to prevent the risk of freezing during transport and storage!</p>
- > Disconnect the >Cooling water inlet< [13] from the water supply.

8.5 Deinstalling the collecting container

PROCEDURE

- Remove the hose from the collecting container.
- > Follow the instructions for the proper disposal of thermal fluid. \rightarrow Page 17, section **»Proper** disposal of resources and consumables«.
- > Dismantle the hose at the >Overflow< [12].



8.6 Uninstalling an external application

PROCEDURE

Disconnect the external application from the temperature control unit.

8.7 Packing

Always use the original packaging! \rightarrow Page 27, section **»Unpacking«**.

8.8 Shipping

NOTE

Temperature control unit transported in a horizontal position

DAMAGE TO THE COMPRESSOR

> Only transport the temperature control unit in an upright position.


NOTE

Temperature control unit transported incorrectly

PROPERTY DAMAGE

- Do not transport by truck on the castors or feet.
- Comply all requirements in this section to avoid damage to the temperature control unit.

Pallet with squared timber for free-standing units

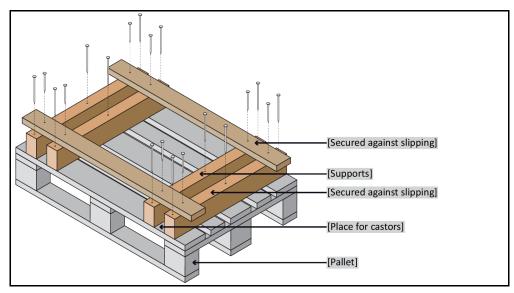

Transport using the lugs, if fitted, on the top of the temperature control unit. Do not transport the temperature control unit alone and without aids.

- Always use the original packaging for transportation.
- Indicate the upright transport position with arrows on the packaging.
- Always transport the temperature control unit upright on a pallet!
- Protect attachments from damage during transportation!
- During transport, place the temperature control unit on squared timber to protect the casters/feet.
- Secure with tensioning belts/lashing straps that are suitable for the weight.
- Additionally secure (depending on model) with plastic film, cardboard and straps.



8.9	Disposal

	Uncontrolled or incorrect opening of the coolant circuit
	RISK OF INJURY AND ENVIRONMENTAL DAMAGE
	 > Work on the coolant circuit and disposal of the refrigerant must be carried out by approved refrigeration/air-conditioning system contractors. > Please strictly observe: → Page 20, section »Temperature control units with fluorinated greenhouse gases/refrigerants«.
ΝΟΤΕ	Improper disposal
	ENVIRONMENTAL DAMAGE
	Spilled/leaked thermal fluid must be discarded immediately and correctly. → Page 17, section »Proper disposal of resources and consumables«.
	To avoid environmental damage, have "disused" temperature control units disposed of exclusively by approved waste management companies (e.g. refrigeration and air conditioning companies).
	Please strictly observe: → Page 20, section »Temperature control units with fluorinated greenhouse gases/refrigerants«.

Huber temperature control units and Huber accessories are made of high quality, recyclable materials. For example: Stainless steel 1.4301 / 1.4401 (V2A), copper, nickel, FKM, Perbunan, NBR, ceramic, carbon, Al-Oxid, red brass, brass, nickel-plated brass and silver solder. Proper recycling of the temperature control unit and accessories can actively help reduce CO₂ emissions in the production of these materials. Follow the laws and regulations of your jurisdiction when disposing material.

8.10 Contact data

INFORMATION

Contact your supplier or local specialist retailer **prior** to returning the temperature control unit. The contact data can be found on our homepage www.huber-online.com under the heading "Contact". Please keep the serial number of the temperature control unit ready. The serial number can be found on the nameplate of the temperature control unit.

8.10.1 Telephone number: Customer Support

If your country is not mentioned in the list below: The responsible service partner can be found on our homepage www.huber-online.com under the heading "Contact".

- Huber Deutschland: +49 781 9603 244
- Huber China: +86 (20) 89001381
- Huber India: +91 80 2364 7966
- Huber Ireland: +44 1773 82 3369
- Huber Italia: +39 0331 181493
- Huber Swiss: +41 (0) 41 854 10 10
- Huber UK: +44 1773 82 3369
- Huber USA: +1 800 726 4877 | +1 919 674 4266
- 8.10.2 Telephone number: Sales

Telephone: +49-781-9603-123

8.10.3 Email address: Customer Support

Email: support@huber-online.com



8.11 Certificate of Compliance

This certificate must be enclosed with the temperature control unit. \rightarrow Page 88, section **»Decontamination/repairs**«.



Annex

Chapter 9

9 Annex

Inspired by **temperature** designed for you

Peter Huber Kältemaschinenbau SE Werner-von-Siemens-Str. 1 77656 Offenburg / Germany

Telefon +49 (0)781 9603-0 Telefax +49 (0)781 57211

info@huber-online.com www.huber-online.com

Technischer Service: +49 (0)781 9603-244

-125 °C ...+425 °C

